



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

Ch
2.22.08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/14/08 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 1/16/08 and 2/14/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

2/15/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 15, day of February, 2008.

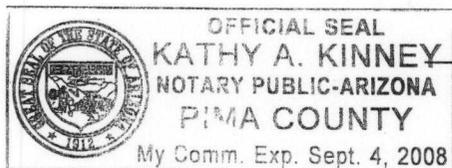

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 75c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

BU
2.22.08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/14/08 TIME 9:50

INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS none

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
02/14/2008

09:52:44

2 of 7 sc
CU
2.22.08

Standard Lot#: 719802I
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:53:33
Diagnostic Test	Pass	09:54:01
Air Blank	0.000	09:54:28
0.100 Cal Check	0.099	09:54:47
Air Blank	0.000	09:55:16
Subject Test	0.000	09:55:46
Air Blank	0.000	09:56:13
Five Minute Wait	ABT*	
Air Blank	0.000	09:56:53

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
02/14/2008 09:58:46

3 of 7sc
CU
2.20-08

Standard Lot#: 719802I
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:59:51
Air Blank	IPS**	10:00:01

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
02/14/2008

10:01:40

4 of 7 SC
lu
2.20.08

Standard Lot#: 719802I
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:02:57
Diagnostic Test	Pass	10:03:25
Air Blank	0.000	10:03:52
0.100 Cal Check	0.100	10:04:11
Air Blank	0.000	10:04:40
Subject Test	0.181*	10:07:49
Air Blank	0.000	10:08:21
Air Blank	0.000	10:08:48
Subject Test	XXX**	10:09:10
Air Blank	0.000	10:09:43

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

02/14/2008

10:10:43

5 of 7 sc
u
2-22-58

Standard Lot#: 719802I

Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:11:10
Cal Check	0.098	10:11:29
Air Blank	0.000	10:11:58

Sarah Cione

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7sc

OH
2-20-08

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 2/14/08 TIME 9:50
Intoxilyzer Serial # 80-005712 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
02/14/2008 09:50:49

7 of 7 SC
ll
2.20-08

Standard Lot#: 719802I
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:51:16
Cal Check	0.099	09:51:35
Air Blank	0.000	09:52:04

Sarah Cione
Operator's Signature