



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
9-25-08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/25/08 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 9/2/08 and 9/25/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


9/25/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 25th, day of September, 2008.

A.M. Hernandez

Notary Public

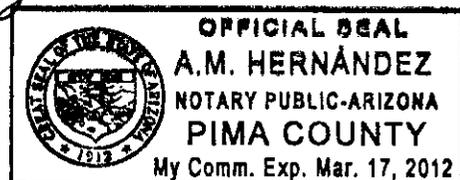

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSY
9.25.08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/25/08 TIME 8:34

INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Sarah Cione

2 of 75c
PH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
09/25/2008 08:36:59

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:37:46
Diagnostic Test	Pass	08:38:13
Air Blank	0.000	08:38:40
0.100 Cal Check	0.098	08:38:59
Air Blank	0.000	08:39:29
Subject Test	0.000	08:39:47
Air Blank	0.000	08:40:14
Five Minute Wait	ABT*	
Air Blank	0.000	08:40:50

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc

104

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
09/25/2008 08:43:00

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:43:42
Diagnostic Test	Pass	08:44:09
Air Blank	0.000	08:44:36
0.100 Cal Check	0.098	08:44:55
Air Blank	0.000	08:45:25
Subject Test	0.127*	08:48:33
Air Blank	0.000	08:49:04
Air Blank	0.000	08:49:31
Subject Test	XXX**	08:50:05
Air Blank	0.000	08:50:37

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

4 of 7 SC

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
09/25/2008 08:51:29

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:51:54
Air Blank	IPS**	08:52:03

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

5 of 7 SC


104

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
09/25/2008 08:52:50

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:53:18
Cal Check	0.000	08:53:37
Air Blank	0.000	08:54:06

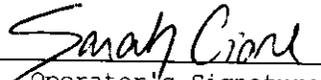

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

B4
9.25.08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/25/08 TIME 8:34
Intoxilyzer Serial # 80-000712 Location DUT HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 SC

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
09/25/2008 08:34:28

BY

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:34:56
Cal Check	0.098	08:35:14
Air Blank	0.000	08:35:44

Sarah Cione
Operator's Signature