



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number PO-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY
11.4.08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/24/08 Cal 1 0.101 AC Cal 2 0.102 AC

Between the dates of 9/25/08 and 10/24/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____
The instrument was returned to service on _____

Sarah Cione

Sarah L. Cione


11/3/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 3rd day of November, 2008.

Art Hernandez

Notary Public.

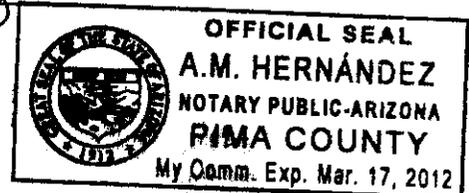

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
11-4-08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/24/08 TIME 17:10

INTOXILYZER SERIAL # 80-000712 LOCATION Vehicle #866

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 102 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7sc

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 90-000712
10/24/2008 17:12:26

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

OAS: S. CIGNE #50616
TUCSON PD

Operator: S. CIGNE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:13:20
Diagnostic Test	Pass	17:13:49
Air Blank	0.000	17:14:16
0.100 Cal Check	0.101	17:14:35
Air Blank	0.000	17:15:04
Subject Test	0.000	17:15:23
Air Blank	0.000	17:15:50
Five Minute Wait	ABT*	
Air Blank	0.000	17:16:25

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7sc

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
10/24/2008 17:16:36

Standard Lot#: 7354016
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test g/210L Time

Air Blank RFI* 17:17:26
Air Blank IPS** 17:17:34

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 750
154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
10/24/2008 17:17:46

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:18:33
Diagnostic Test	Pass	17:19:02
Air Blank	0.000	17:19:29
0.100 Cal Check	0.101	17:19:48
Air Blank	0.000	17:20:17
Subject Test	0.053*	17:23:26
Air Blank	0.000	17:23:58
Air Blank	0.000	17:24:25
Subject Test	XXX**	17:24:44
Air Blank	0.000	17:25:17

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test-Sequence

5 of 7sc
184

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
10/24/2008 17:25:45

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616 #
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:26:12
Cal Check	0.102	17:26:31
Air Blank	0.000	17:27:01

Sarah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

154
11.4.08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/24/08 TIME 17:10
Intoxilyzer Serial # 80-000712 Location Vehicle # 866

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 75c

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
10/24/2008 17:10:44

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:11:12
Cal Check	0.101	17:11:30
Air Blank	0.000	17:12:00

Sarah Cione

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2 sc

134
11-4-08

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 10/05/08
 QAS: CLEARY +44524
 Intoxilyzer Serial Number: 80-000712
 New Lot Number: 7354101G
 New Expiration Date: 12/28/09

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press
 ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE San Cleary +44524

2 of 2 sc

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
10/05/2008 18:46:34

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

DAS: S. CLEARY #44524
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	18:47:01
Cal Check	0.101	18:47:20
Air Blank	0.000	18:47:50


Operator's Signature