



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
2-17-09

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/12/09 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 1/18/09 and 2/12/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
 The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


2/17/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17th, day of February, 20 09.

A.M. Hernández

Notary Public

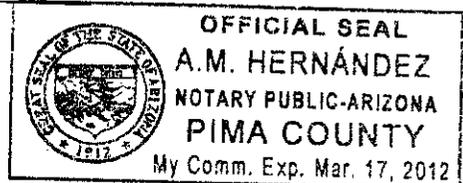

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
2-17-09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/12/09 TIME 7:17

INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
02/12/2009

07:18:58

2 of 7 SC

TS4

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:19:43
Diagnostic Test	Pass	07:20:11
Air Blank	0.000	07:20:38
0.100 Cal Check	0.100	07:20:57
Air Blank	0.000	07:21:26
Subject Test	0.000	07:21:43
Air Blank	0.000	07:22:11
Five Minute Wait	ABT*	
Air Blank	0.000	07:22:45

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
02/12/2009

07:23:30

3 of 7 SC
184

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:23:55
Air Blank	IPS**	07:24:04

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
02/12/2009

07:24:48

4 of 7 SC
124

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:25:36
Diagnostic Test	Pass	07:26:03
Air Blank	0.000	07:26:31
0.100 Cal Check	0.099	07:26:49
Air Blank	0.000	07:27:19
Subject Test	0.062*	07:30:27
Air Blank	0.000	07:30:58
Air Blank	0.000	07:31:26
Subject Test	XXX**	07:31:56
Air Blank	0.000	07:32:29

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
02/12/2009

07:33:13

5 of 7 SC

134

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50610 -
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:33:41
Cal Check	0.099	07:33:59
Air Blank	0.000	07:34:29

Sarah Cione

Operator's signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
2-17-09

QA Specialist Sarah Cone AGENCY Tucson Police Department
DATE 2/12/09 TIME 7:17
Intoxilyzer Serial # 80-000712 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cone

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
02/12/2009

07:17:29

7 of 7 sc

Standard Lot#: 735401G
Last Changed By: S. CLEARY #44524

154

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:17:56
Cal Check	0.100	07:18:15
Air Blank	0.000	07:18:44

Sarah Cione
Operator's Signature