



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*By
8-31-2009*

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 08/27/09 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 07/28/09 and 08/27/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


08/31/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 31, day of August, 2009.

A.M. Hernandez

Notary Public

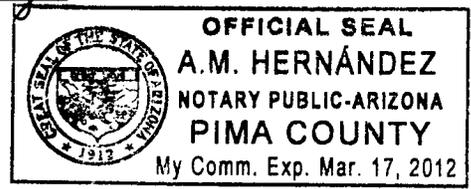

*Includes documents specified in A.A.C. R13-10-104.A.5

187 CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

BY
8.31.09

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl Steaub AGENCY Tucson Police Department

DATE 08/27/09 TIME 08:26

INTOXILYZER SERIAL # 80-000712 LOCATION CRIME LAB

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Steaub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

287
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
08/27/2009

08:29:14

Standard Lot#: 820501I
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:30:08
Diagnostic Test	Pass	08:30:36
Air Blank	0.000	08:31:03
0.100 Cal Check	0.099	08:31:22
Air Blank	0.000	08:31:51
Subject Test	0.000	08:32:16
Air Blank	ABT*	08:32:23
Air Blank	0.000	08:32:51

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
08/27/2009

08:33:00

Standard Lot#: 820501I
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:33:36
Air Blank	IPS**	08:33:44

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

487
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
08/27/2009

08:34:05

PM

Standard Lot#: 820501I
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:35:04
Diagnostic Test	Pass	08:35:32
Air Blank	0.000	08:35:59
0.100 Cal Check	0.099	08:36:18
Air Blank	0.000	08:36:47
Subject Test	0.217*	08:39:56
Air Blank	0.000	08:40:28
Air Blank	0.000	08:40:55
Subject Test	XXX**	08:41:23
Air Blank	0.000	08:41:57

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

507
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
08/27/2009 08:42:19

104

Standard Lot#: 820501I
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:42:46
Cal Check	0.099	08:43:05
Air Blank	0.000	08:43:34



Operator's Signature

6 of 7
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

BH
8.31.09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl Stramb AGENCY Tucson Police Department
DATE 08/27/09 TIME 08:26
Intoxilyzer Serial # 80-000712 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 0.99 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Stramb

707
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

08/27/2009

08:26:08

Standard Lot#: 820501I


Last Changed By: C. STRAUB #50772 .

QAS: C. STRAUB #50772

TUCSON PD

SM

Test	g/210L	Time
Air Blank	0.000	08:26:36
Cal Check	0.099	08:26:54
Air Blank	0.000	08:27:24

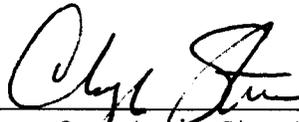

Operator's Signature