



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

By
3-10-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/4/2010 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 2/4/2010 and 3/4/2010 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


3/4/2010

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 4th, day of March, 2010.

A.M. Hernández

Notary Public

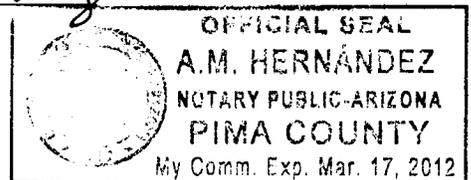

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 9 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

Tsh
3-10-2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 3/4/2010 TIME 1322
INTOXILYZER SERIAL # 80-000712 LOCATION Vehicle #437

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 9 sc
134

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
03/04/2010 13:24:17

Standard Lot#: 8205011
Last Changed By: C. STRAUB #50772

DAS: S. CLONE #50616
TUCSON PD

Operator: S. CLONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:25:10
Diagnostic Test	Pass	13:25:39
Air Blank	0.000	13:26:06
0.100 Cal Check	0.098	13:26:25
Air Blank	0.000	13:26:54
Subject Test	0.000	13:27:26
Air Blank	0.000	13:27:53
Five Minute Wait	ABT*	
Air Blank	0.000	13:28:34

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 9 sc

124

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
03/04/2010 13:28:57

Standard Lot#: 8205011
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	13:29:35
Air Blank	IPS**	13:29:45

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 9 sc

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
03/04/2010 13:29:54

Standard Lot#: 8205011
Last Changed By: C. STRAUB #50772

DAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:30:45
Diagnostic Test	Pass	13:31:13
Air Blank	0.000	13:31:41
0.100 Cal Check	0.098	13:31:59
Air Blank	0.000	13:32:29
Subject Test	0.109*	13:35:37
Air Blank	0.000	13:36:09
Air Blank	0.000	13:36:37
Subject Test	XXX**	13:36:56
Air Blank	0.000	13:37:29

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 9 sc

124

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
03/04/2010 13:37:49

Standard Lot#: 8205011
Last Changed By: C. STRAUB #50772

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:38:16
Cal Check	0.099	13:38:35
Air Blank	0.000	13:39:05

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 9 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
3.10-2010

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/4/2010 TIME 1322
Intoxilyzer Serial # 80-000712 Location Vehicle # 437

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 9 sc
154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
03/04/2010 13:22:42

Standard Lot#: 8205011
Last Changed By: C. STRAUB #50772

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:23:10
Cal Check	0.099	13:23:29
Air Blank	0.000	13:23:58

Sarah Cione
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 sc

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TSC
3.10.2010

Date: 3/4/2010

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000712

New Lot Number: 642845

New Expiration Date: 7/23/2011

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
 - Enter Target Value: 0.100
 - Enter Control Lot #: From the tank label
 - Enter Exp Date: From the tank label (mm/dd/yyyy)
 - Edit info Y/N? Y or N
 - Update Std Info? Y ENTER
 - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

9 of 9 sc

PH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
03/04/2010 13:43:49

Standard Lot#: 642845
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:44:16
Cal Check	0.100	13:44:35
Air Blank	0.000	13:45:05

Sarah Cione
Operator's Signature