



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

Handwritten note: B4 8.23.2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/19/10 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 7/22/10 and 8/19/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Signature of Sarah L. Cione

Sarah L. Cione


8/23/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 23 day of August, 2010.

Signature of Notary Public A.M. Hernandez

Notary Public

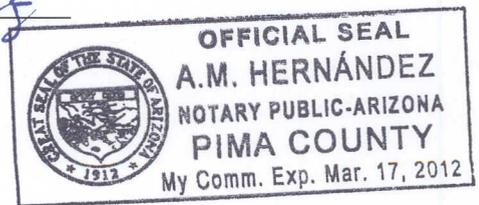

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
8.23.2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 8/19/10 TIME 942

INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
08/19/2010

09:45:41

2 of 7sc

BM

Standard Lot#: 719919
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:46:37
Diagnostic Test	Pass	09:47:05
Air Blank	0.000	09:47:32
0.100 Cal Check	0.100	09:47:51
Air Blank	0.000	09:48:20
Subject Test	0.000	09:49:01
Air Blank	0.000	09:49:28
Five Minute Wait	ABT*	
Air Blank	0.000	09:50:07

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
08/19/2010

09:50:25

3 of 7_{sc}

PH

Standard Lot#: 719919
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:50:57
Air Blank	IPS**	09:51:06

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
08/19/2010 09:52:59

SM

Standard Lot#: 719919
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:54:05
Diagnostic Test	Pass	09:54:33
Air Blank	0.000	09:55:00
0.100 Cal Check	0.098	09:55:19
Air Blank	0.000	09:55:48
Subject Test	0.043*	09:58:57
Air Blank	0.000	09:59:28
Air Blank	0.000	09:59:55
Subject Test	XXX**	10:00:20
Air Blank	0.000	10:00:53

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
08/19/2010 10:01:31

5 of 7_{sc}

SM

Standard Lot#: 719919
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:01:58
Cal Check	0.098	10:02:17
Air Blank	0.000	10:02:46

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

BY
8.23.2010

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 8/19/10 TIME 942
Intoxilyzer Serial # 80-000712 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
08/19/2010 09:42:49

7 of 7 sc
BY

Standard Lot#: 719919
Last Changed By: S. CLEARY #44524

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:43:16
Cal Check	0.100	09:43:35
Air Blank	0.000	09:44:04

Sarah Cione
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2_{sc}

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

By
8.23.2010

Date: 08/01/10
 QAS: SEAN CLEARY
 Intoxilyzer Serial Number: 80-000712
 New Lot Number: 719919
 New Expiration Date: 02/02/12

- (X) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (X) 2. Remove the old gas tank and install the new one.
- (X) 3. Hit the "ESC" key twice and swipe the QAS card.
- (X) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (X) 5. Hit the "ESC" key twice.
- (X) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press ENTER.
- (X) 7. Air blank completed.
- (X) 8. Calibration check completed. Test results 0. 099 AC.
- (X) 9. Air blank completed.
- (X) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sean Cleary #445201

2 of 2 sc


154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000712
08/01/2010 19:28:42

Standard Lot#: 719919
Last Changed By: S. CLEARY #44524

QAS: S. CLEARY #44524
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	19:29:09
Cal Check	0.099	19:29:28
Air Blank	0.000	19:29:57

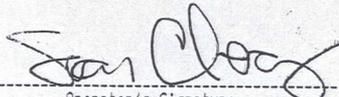
Operator's Signature