



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

Tsh  
5-2-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-28-11 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 3-31-11 and 4-28-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


4.28-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 28, day of April, 20 11.

Isabel C. Trujillo

Notary Public

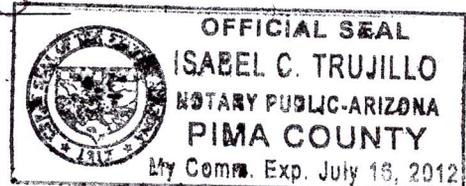

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

1078  
on  
154  
5.2.2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department  
DATE 4-28-11 TIME 0931  
INTOXILYZER SERIAL # 80-000712 LOCATION ODW DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley  
\_\_\_\_\_

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
04/28/2011

09:32:37

2008  
CM  
104

Standard Lot#: 719919  
Last Changed By: S. CLEARY #44524

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:33:39
Diagnostic Test	Pass	09:34:07
Air Blank	0.000	09:34:34
0.100 Cal Check	0.098	09:34:53
Air Blank	0.000	09:35:23
Subject Test	RFI*	09:36:24
Air Blank	IPS**	09:36:41

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
04/28/2011 09:36:57

3-1-08  
CH  
TS4

Standard Lot#: 719919  
Last Changed By: S. CLEARY #44524

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:37:54
Diagnostic Test	Pass	09:38:22
Air Blank	0.000	09:38:49
0.100 Cal Check	0.099	09:39:08
Air Blank	0.000	09:39:37
Subject Test	0.000	09:39:53
Air Blank	ABT*	09:40:00
Air Blank	0.000	09:40:28

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
04/28/2011 09:40:58

4.168  
cu  
B4

Standard Lot#: 719919  
Last Changed By: S. CLEARY #44524

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:41:49
Diagnostic Test	Pass	09:42:16
Air Blank	0.000	09:42:44
0.100 Cal Check	0.099	09:43:02
Air Blank	0.000	09:43:32
Subject Test	0.095*	09:46:40
Air Blank	0.000	09:47:12
Air Blank	0.000	09:47:40
Subject Test	XXX**	09:48:06
Air Blank	0.000	09:48:39

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
04/28/2011 09:49:10

508  
CM

-----  
Standard Lot#: 719919  
Last Changed By: S. CLEARY #44524  
-----

BY

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:49:38
Cal Check	0.098	09:49:56
Air Blank	0.000	09:50:26

Claire Conley  
Operator's Signature

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
04/28/2011 09:51:42

*b. b. g.  
cu*

-----  
Standard Lot#: 719919  
Last Changed By: S. CLEARY #44524  
-----

*AM*

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:52:09
Cal Check	0.099	09:52:28
Air Blank	0.000	09:52:57

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

*accidentally ran*

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7088  
CM  
134  
5.2.2011

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 4-28-11 TIME 0931  
Intoxilyzer Serial # 80-000712 Location ODW DM HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
04/28/2011 09:31:07

8/28  
CW

-----  
Standard Lot#: 719919  
Last Changed By: S. CLEARY #44524  
-----

BY

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:31:35
Cal Check	0.098	09:31:53
Air Blank	0.000	09:32:23

*Clare Conley*

\_\_\_\_\_  
Operator's Signature