



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TSh
 7-24-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/19/12 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 6/21/12 and 7/19/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


7/20/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 20th, day of July, 2012.

Elma Escarcega

Notary Public

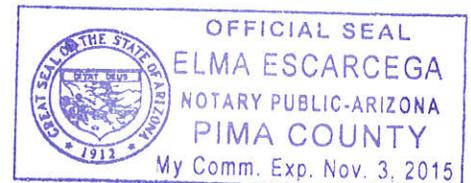

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 SC
TS4
7.24.2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/19/12 TIME 1125

INTOXILYZER SERIAL # 80-000712 LOCATION DUIHQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
07/19/2012

11:27:50

2 of 7_{sc}
BSH

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:28:31
Diagnostic Test	Pass	11:29:05
Air Blank	0.000	11:29:32
0.100 Cal Check	0.098	11:29:53
Air Blank	0.000	11:30:23
Subject Test	0.000	11:30:43
Air Blank	0.000	11:31:10
Wait	ABT*	
Air Blank	0.000	11:32:27

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
07/19/2012 11:33:33

3 of 7_{sc}
B4

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:34:00
Air Blank	IPS**	11:34:07

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
07/19/2012

11:38:06

4 of 7 sc

184

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:38:54
Diagnostic Test	Pass	11:39:28
Air Blank	0.000	11:39:55
0.100 Cal Check	0.097	11:40:16
Air Blank	0.000	11:40:46
Subject Test	0.202*	11:43:53
Air Blank	0.000	11:44:27
Air Blank	0.000	11:44:54
Subject Test	XXX**	11:45:23
Air Blank	0.000	11:45:57

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
07/19/2012 11:48:13

5 of 7 sc
JSh

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:48:45
0.100 Cal Check	0.098	11:49:06
Air Blank	0.000	11:49:36

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

6 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
7.24.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 7/19/12 TIME 1125
Intoxilyzer Serial # 80-000712 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
07/19/2012 11:25:29

7 of 7 sc
TSH

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:25:58
0.100 Cal Check	0.098	11:26:18
Air Blank	0.000	11:26:48

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2 sc

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TS4
7-24-2012

Date: 6/23/12
 QAS: J HAWKINS 42401
 Intoxilyzer Serial Number: 80-000712
 New Lot Number: 2351100A2
 New Expiration Date: 11/1/2013

1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
2. Remove the old gas tank and install the new one.
3. Hit the "ESC" key twice and swipe the QAS card.
4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
5. Hit the "ESC" key twice.
6. Go to the "M-Maintenance" menu and select "Control Testing." Select "D" for dry control test. After selection is made press ENTER.
7. Air blank completed.
8. Calibration check completed. Test results 0. 099 AC.
9. Air blank completed.
10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE [Signature]

2 of 2 sc

154

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
06/23/2012 01:02:35

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CLONE #50616
TUCSON PD

Operator: J. HAWKINS #42401
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	01:03:07
0.100 Cal Check	0.099	01:03:28
Air Blank	0.000	01:03:57

Preventative Maintenance NOT
Performed -- Timer NOT Reset

