



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

154
11-4-08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/23/08 Cal 1 0.102 AC Cal 2 0.103 AC

Between the dates of 9/29/08 and 10/23/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


11/3/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 3rd, day of November, 2008.

A.M. Hernandez

Notary Public

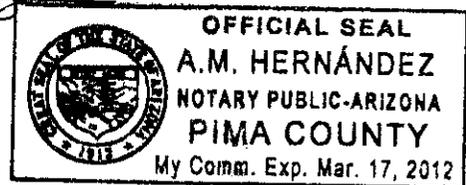

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

134
11.4.08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/23/08 TIME 8:38

INTOXILYZER SERIAL # 80-001108 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 103 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 8_{sc}

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
10/23/2008

08:40:58

154

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:41:53
Diagnostic Test	Pass	08:42:21
Air Blank	0.000	08:42:48
0.100 Cal Check	0.103	08:43:07
Air Blank	0.000	08:43:37
Subject Test	0.000	08:44:19
Air Blank	0.000	08:44:46
Five Minute Wait	ABT*	
Air Blank	0.000	08:46:10

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8
SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
10/23/2008 08:46:32

154

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:47:28
Air Blank	0.000	08:47:56

*RFI Detect

Not a Successfully Completed Test Sequence

4 of 8 sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
10/23/2008

08:48:59

154

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	08:49:33
Air Blank	0.000	08:50:02

*Improper Sample

Not a Successfully Completed Test Sequence

5 of 8
sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
10/23/2008

08:51:56

154

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:52:49
Diagnostic Test	Pass	08:53:16
Air Blank	0.000	08:53:44
0.100 Cal Check	0.103	08:54:02
Air Blank	0.000	08:54:32
Subject Test	0.153*	08:57:41
Air Blank	0.000	08:58:12
Air Blank	0.000	08:58:40
Subject Test	XXX**	08:59:00
Air Blank	0.000	08:59:33

*Deficient Sample
 **Invalid Sample
 Begin new deprivation period

Not a Successfully Completed Test Sequence

6 of 8 5c

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
10/23/2008 09:00:47

154

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:01:14
Cal Check	0.103	09:01:33
Air Blank	0.000	09:02:02

Sarah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

154
11-4-08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/23/08 TIME 8:38
Intoxilyzer Serial # 80-001108 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

8 of 8
SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
10/23/2008 08:38:48

TS4

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:39:15
Cal Check	0.102	08:39:34
Air Blank	0.000	08:40:03

Sarah Cione
Operator's Signature