



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY
12-22-08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/18/08 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 11/20/08 and 12/18/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


12/22/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22, day of December, 2008.

A.M. Hernandez

Notary Public

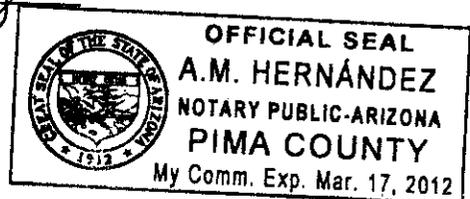

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TS4
12-22-08

QA SPECIALIST Sarah Ciore AGENCY Tucson Police Department
DATE 12/18/08 TIME 9:32
INTOXILYZER SERIAL # 80-001108 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Ciore

2 of 7_{SC}

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
12/18/2008

09:34:14

154

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:35:01
Diagnostic Test	Pass	09:35:29
Air Blank	0.000	09:35:56
0.100 Cal Check	0.102	09:36:15
Air Blank	0.000	09:36:44
Subject Test	0.000	09:37:15
Air Blank	0.000	09:37:43
Five Minute Wait	ABT*	
Air Blank	0.000	09:38:17

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
12/18/2008

09:38:23

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

TS4

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:38:48
Air Blank	IPS**	09:38:57

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 75C

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
12/18/2008

09:39:45

B4

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:40:54
Diagnostic Test	Pass	09:41:21
Air Blank	0.000	09:41:48
0.100 Cal Check	0.102	09:42:07
Air Blank	0.000	09:42:37
Subject Test	0.178*	09:45:45
Air Blank	0.000	09:46:17
Air Blank	0.000	09:46:45
Subject Test	XXX**	09:47:12
Air Blank	0.000	09:47:45

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 75c

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
12/18/2008 09:48:00

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

754

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:48:27
Cal Check	0.100	09:48:46
Air Blank	0.000	09:49:15

Sarah Cione

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 75c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TSC
12-22-08

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 12/18/08 TIME 9:32
Intoxilyzer Serial # 80-001108 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
12/18/2008

09:32:35

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
- TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:33:02
Cal Check	0.101	09:33:21
Air Blank	0.000	09:33:51

Sarah Cione

Operator's Signature

7 of 7sc

134