



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
4-7-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/31/11 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 3/3/11 and 3/31/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


4/6/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 6<sup>th</sup> day of April, 20 11.

A.M. Hernandez

Notary Public

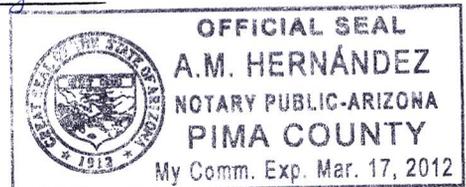

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 7<sub>sc</sub>

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

124  
4.7.2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/31/11 TIME 1443

INTOXILYZER SERIAL # 80-001108 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7sc

134

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
03/31/2011 14:45:32

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:46:19
Diagnostic Test	Pass	14:46:47
Air Blank	0.000	14:47:14
0.100 Cal Check	0.101	14:47:33
Air Blank	0.000	14:48:02
Subject Test	0.000	14:48:21
Air Blank	0.000	14:48:49
Five Minute Wait	ABT*	
Air Blank	0.000	14:49:59

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7sc

1704

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
03/31/2011 14:50:08

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:50:39
Air Blank	IPS**	14:50:47

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7 SC

TS4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
03/31/2011 14:50:52

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:51:38
Diagnostic Test	Pass	14:52:06
Air Blank	0.000	14:52:33
0.100 Cal Check	0.103	14:52:52
Air Blank	0.000	14:53:21
Subject Test	0.025*	14:56:29
Air Blank	0.000	14:57:01
Air Blank	0.000	14:57:28
Subject Test	XXX**	14:57:43
Air Blank	0.000	14:58:16

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7 sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
03/31/2011 14:59:05

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

154

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:59:33
Cal Check	0.101	14:59:51
Air Blank	0.000	15:00:21

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

6 of 75c

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSY  
4.7.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/31/11 TIME 1443  
Intoxilyzer Serial # 80-001108 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
03/31/2011 14:43:14

7 of 7 sc

AB4

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:43:41
Cal Check	0.101	14:44:00
Air Blank	0.000	14:44:29

*Sarah Cione*

Operator's Signature