



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
9.28.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/18/11 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 7/22/11 and 8/18/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


9/22/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22, day of September, 20 11.

A.M. Hernandez

Notary Public

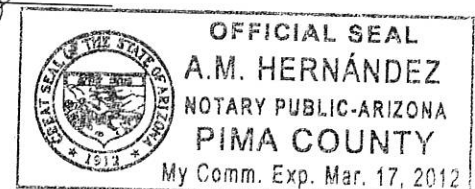

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

BS4
9.28.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 8/18/11 TIME 1533
INTOXILYZER SERIAL # 80-001108 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
08/18/2011

15:35:24

2 of 8_{SC}
BM

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:36:24
Diagnostic Test	Pass	15:36:51
Air Blank	0.000	15:37:18
0.100 Cal Check	0.099	15:37:37
Air Blank	0.000	15:38:07
Subject Test	0.000	15:38:52
Air Blank	0.000	15:39:20
Five Minute Wait	ABT*	
Air Blank	0.000	15:39:57

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
08/18/2011

15:41:40

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	15:42:33
Air Blank	IPS**	15:42:44

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
08/18/2011

15:43:16

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:44:18
Diagnostic Test	Pass	15:44:45
Air Blank	0.000	15:45:12
0.100 Cal Check	0.100	15:45:31
Air Blank	0.000	15:46:01
Subject Test	0.087*	15:49:09
Air Blank	0.000	15:49:40
Air Blank	0.000	15:50:08
Subject Test	0.123**	15:50:39
Air Blank	0.000	15:51:16
Five Minute Wait		
Air Blank	0.000	15:56:07
Subject Test	0.000**	15:56:47
Air Blank	0.000	15:57:15

*Deficient Sample

**No .020 agreement

Mouth alcohol sample analyzed as an actual sample. Waited 5 min. to demonstrate evaporation.

SC

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
08/18/2011

15:57:53

5 of 8 SC

by

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:58:42
Diagnostic Test	Pass	15:59:10
Air Blank	0.000	15:59:37
0.100 Cal Check	0.099	15:59:56
Air Blank	0.000	16:00:25
Subject Test	XXX*	16:00:56
Air Blank	0.000	16:01:29

*Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
08/18/2011

16:03:21

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:03:48
Cal Check	0.099	16:04:07
Air Blank	0.000	16:04:37



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 8/18/11 TIME 1533
Intoxilyzer Serial # 80-001108 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
08/18/2011

15:33:24

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:33:52
Cal Check	0.101	15:34:11
Air Blank	0.000	15:34:40

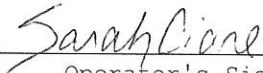

Operator's Signature

8 of 8 sc

154