



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

je

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-28-2006 Cal 1 0.098 AC Cal 2 0.096 AC

Between the dates of 12-1-2006 and 12-28-2006 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .


Terry S. Gallegos

Terry S. Gallegos

1.4.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.

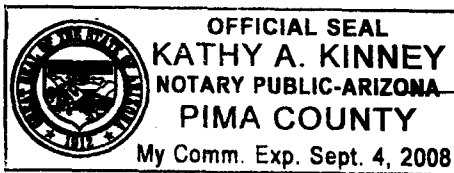

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

10 of 10
Pg

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 20847 AGENCY Tucson Police Department

DATE 12-28-2006 TIME 1756

INTOXILYZER SERIAL # 80-001163 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning. ✓✓✓
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001163
12/28/2006

17:58:29

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:59:24
Diagnostic Test	Pass	17:59:51
Air Blank	0.000	18:00:19
0.100 Cal Check	0.097	18:00:38
Air Blank	0.000	18:01:07
Subject Test	0.000	18:01:24
Air Blank	0.000	18:01:51
Five Minute Wait	ABT*	
Air Blank	0.000	18:02:56

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001163
12/28/2006

18:03:06

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:03:58
Diagnostic Test	Fail*	18:04:02
Air Blank	0.000	18:04:34

* Voltage/Current Test

RFI during diagnostic test B4 28847

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001163

12/28/2006

18:05:53

Standard Lot#: 502602I

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:06:41
Diagnostic Test	Pass	18:07:08
Air Blank	0.000	18:07:36
0.100 Cal Check	0.097	18:07:54
Air Blank	0.000	18:08:24
Subject Test	RFI*	18:08:44
Air Blank	0.000	18:09:13

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001163
12/28/2006

18:09:33

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	18:10:11
Air Blank	0.000	18:10:40

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001163
12/28/2006

18:12:36

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:13:25
Diagnostic Test	Pass	18:13:52
Air Blank	0.000	18:14:20
0.100 Cal Check	0.098	18:14:39
Air Blank	0.000	18:15:08
Subject Test	0.093*	18:18:16
Air Blank	0.000	18:18:48
Air Blank	0.000	18:19:15
Subject Test	DEF**	18:22:24
Air Blank	0.000	18:22:52

*Deficient Sample

**Deficient Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001163

12/28/2006

18:24:31

Standard Lot#: 502602I

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:25:29
Diagnostic Test	Pass	18:25:56
Air Blank	0.000	18:26:23
0.100 Cal Check	0.098	18:26:42
Air Blank	0.000	18:27:11
Subject Test	XXX*	18:27:37
Air Blank	0.000	18:28:11

*Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001163

12/28/2006

18:28:39

8 of 10
B4

Standard Lot#: 502602I

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	18:29:07
Cal Check	0.096	18:29:25
Air Blank	0.000	18:29:55


Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

9 of 10
184

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 12-29-2006 TIME 1756
Intoxilyzer Serial # 80-001163 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001163

12/28/2006

17:56:53

100810
P4

Standard Lot#: 502602I

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:57:21
Cal Check	0.098	17:57:39
Air Blank	0.000	17:58:09


Operator's Signature