



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

WMB
5-2581

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4.26.2007 Cal 1 0.099 AC Cal 2 0.100 AC

Between the dates of 3.26.2007 and 4.26.2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry S. Gallegos

5.24.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24, day of May, 2007.

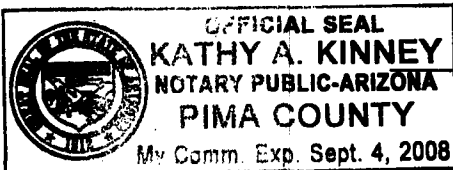

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police DepartmentDATE 4-26-2007 TIME 1835INTOXILYZER SERIAL # 80-001163 LOCATION DUI HQ(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (
- ☒
-) 1. Clock time check.
-
- (
- ☒
-) 2. Date check.

OPERATIONAL TESTS

- (
- ☒
-) 1. Alcohol-free subject test result 0.
- 000
- AC.
-
- (
- ☒
-) 2. Error recognition logic system functioning.
-
- Not a Successfully Completed Test Sequence printed ✓✓✓
-
- (
- ☒
-) 3. Proper sample recognition system.
-
- Not a Successfully Completed Test Sequence printed
-
- Deficient sample printed.
-
- (
- ☒
-) 4. Standard Calibration Check standard 0.
- 100
- AC
-
- Result: 0.
- 100
- AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

1 of 7
TS4 7/15
5-2507

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001163

04/26/2007

18:37:24

2087
TS4 5-2507

Standard Lot#: 635301I

Last Changed By: G. EPPLEY

#26227

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:38:12
Diagnostic Test	Pass	18:38:39
Air Blank	0.000	18:39:07
0.100 Cal Check	0.100	18:39:26
Air Blank	0.000	18:39:55
Subject Test	0.009	18:40:12
Air Blank	ABT*	18:40:26
Air Blank	0.000	18:40:55

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001163

04/26/2007

18:41:07

3 of 7
TS4 5-2507

Standard Lot#: 635301I

Last Changed By: G. EPPLEY

#26227

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:41:59
Diagnostic Test	Pass	18:42:27
Air Blank	0.000	18:42:54
0.100 Cal Check	0.100	18:43:13
Air Blank	0.000	18:43:42
Subject Test	0.000	18:44:00
Air Blank	RFI*	18:44:10
Air Blank	IPS**	18:44:18

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001163

04/26/2007

18:44:28

Standard Lot#: 635301I

Last Changed By: G. EPPLEY

#26227

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:45:18
Diagnostic Test	Pass	18:45:45
Air Blank	0.000	18:46:13
0.100 Cal Check	0.099	18:46:31
Air Blank	0.000	18:47:01
Subject Test	0.074*	18:50:09
Air Blank	0.000	18:50:42
Air Blank	0.000	18:51:09
Subject Test	XXX**	18:51:28
Air Blank	0.000	18:52:02

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4087
B4 5-28-07

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001163

04/26/2007

18:53:23


Standard Lot#: 635301I

Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	18:53:50
Cal Check	0.100	18:54:09
Air Blank	0.000	18:54:38


Operator's Signature

50f7
Ry 5/2/07

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

6077
B4 5/2007

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 4-26-2007 TIME 1835
Intoxilyzer Serial # 80-001163 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE



DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001163

04/26/2007

18:35:50

Standard Lot#: 635301I

Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	18:36:17
Cal Check	0.099	18:36:36
Air Blank	0.000	18:37:05


Operator's Signature

7af7
B4^{me}₅₂₅₀₇