



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-00116⁵³ maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY
3/17/08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/13/08 Cal 1 0.101 AC Cal 20.100 AC

Between the dates of 2/8/08 and 3/13/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

3/14/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14, day of March, 2008.

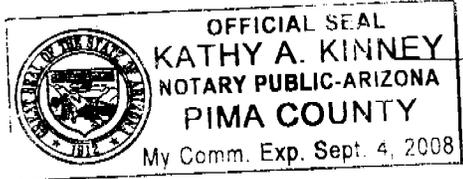

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 85C
TS4
3.17.08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/13/08 TIME 10:08

INTOXILYZER SERIAL # 80-001163 LOCATION Team 2 - DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS G2 was performed greater than 31 days.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
03/13/2008 10:10:18

2 of 8 SC
154

Standard Lot#: 732503I
Last Changed By: T. O HARA #36676

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:11:06
Diagnostic Test	Pass	10:11:33
Air Blank	0.000	10:12:00
0.100 Cal Check	0.101	10:12:19
Air Blank	0.000	10:12:48
Subject Test	0.000	10:13:04
Air Blank	0.000	10:13:31
Five Minute Wait	ABT*	
Air Blank	0.000	10:14:13

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
03/13/2008 10:14:18

3 of 85c
TSH

Standard Lot#: 732503I
Last Changed By: T. O HARA #38876

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:15:04
Diagnostic Test	Pass	10:15:32
Air Blank	RFI*	10:15:47
Air Blank	0.000	10:16:15

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
03/13/2008 10:17:37

4 of 8 sc
TS4

Standard Lot#: 732503I
Last Changed By: T. O HARA #38876

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	10:18:29
Air Blank	0.000	10:18:57

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
03/13/2008 10:19:02

5 of 8 sc
BM

Standard Lot#: 732503I
Last Changed By: T. O HARA #38876

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:19:58
Diagnostic Test	Pass	10:20:25
Air Blank	0.000	10:20:53
0.100 Cal Check	0.101	10:21:12
Air Blank	0.000	10:21:41
Subject Test	0.135*	10:24:49
Air Blank	0.000	10:25:21
Air Blank	0.000	10:25:49
Subject Test	XXX**	10:26:14
Air Blank	0.000	10:26:48

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
03/13/2008 10:27:21

6 of 8 sc
SA

Standard Lot#: 732503I
Last Changed By: T. O HARA #38876

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:27:49
Cal Check	0.100	10:28:07
Air Blank	0.000	10:28:37

Sarah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8 sc
154

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/13/08 TIME 10:08
Intoxilyzer Serial # 80-001163 Location Team 2 - DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
03/13/2008 10:08:26

8 of 85c
134

Standard Lot#: 732503T
Last Changed By: T. O HARA #38876

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:08:53
Cal Check	0.101	10:09:12
Air Blank	0.000	10:09:41

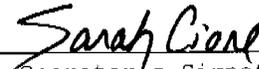

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2
SC

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TS4
3-17-08

Date: 3-6-08
 QAS: OHARA
 Intoxilyzer Serial Number: 80-001163
 New Lot Number: 732503I
 New Expiration Date: 26NOV2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 102 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE OHARA #38876

2 of 2 sc

TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
03/06/2008 19:39:04

Standard Lot#: 7325031
Last Changed By: T. O HARA #38876

QAS: T. O HARA #38876
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	19:39:31
Cal Check	0.102	19:39:50
Air Blank	0.000	19:40:19

T. O Hara # 38876
Operator's Signature
