



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

SC
4/27/10

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-201163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1-15-09 Cal 1 0.100 AC Cal 2 0.097 AC

Between the dates of 12-12-08 and 1-15-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


4-27-10

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 27, day of April, 2010.

Sue A. Soza

Notary Public

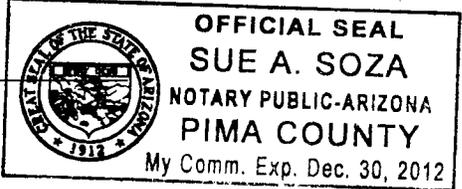

*Includes documents specified in A.A.C. R13-10-104.A.5



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

134
2-4-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1-15-09 Cal 1 0.100 AC Cal 2 0.097 AC

Between the dates of 12-12-08 and 1-15-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


2-3-09

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 3rd day of February, 2009.

A. M. Hernandez

Notary Public

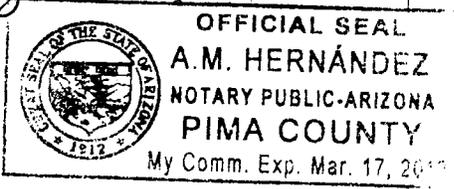

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 25c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
2.4.09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 12/29/08 TIME 15:05
Intoxilyzer Serial # 80-001163 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

2 of 25c

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
12/29/2008 15:05:44

TS4

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:06:12
Cal Check	0.099	15:06:30
Air Blank	0.000	15:07:00

Sarah Cione
Operator's Signature

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1587
cu
124
2-4-09

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 1-15-09 TIME 0910

INTOXILYZER SERIAL # 80-001163 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
01/15/2009

09:12:16

267
on
BH

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:13:02
Diagnostic Test	Pass	09:13:30
Air Blank	0.000	09:13:57
0.100 Cal Check	0.099	09:14:16
Air Blank	0.000	09:14:45
Subject Test	0.000	09:15:06
Air Blank	ABT*	09:15:13
Air Blank	0.000	09:15:41

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
01/15/2009

09:15:48

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:16:52
Air Blank	IPS**	09:17:02

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

3-17
en

TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
01/15/2009

09:17:51

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:18:42
Diagnostic Test	Pass	09:19:10
Air Blank	0.000	09:19:37
0.100 Cal Check	0.099	09:19:56
Air Blank	0.000	09:20:25
Subject Test	0.361*	09:23:34
Air Blank	0.000	09:24:08
Air Blank	0.000	09:24:35
Subject Test	XXX**	09:25:09
Air Blank	0.000	09:25:42

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
01/15/2009 09:26:06

567
en

154

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:26:34
Cal Check	0.097	09:26:52
Air Blank	0.000	09:27:22

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6-57
aw
734
2/4/09

QA Specialist Clare Conley AGENCY Tucson Police Department

DATE 1-15-09 TIME 0910

Intoxilyzer Serial # 80-001163 Location ODW DM HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-001163

01/15/2009

09:10:55

Standard Lot#: 735401G

Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517

• TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:11:22
Cal Check	0.100	09:11:41
Air Blank	0.000	09:12:11

Clare Conley

Operator's Signature

727
cu

TS4