



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 82-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*TSY
2-12-09*

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-12-09 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 1-15-09 and 2-12-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


2-12-09

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 12th, day of February, 2009.

A.M. Hernandez

Notary Public

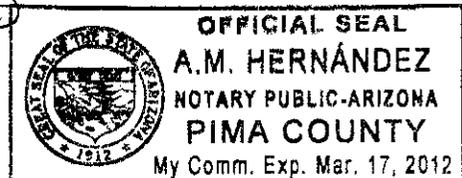

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1009
CR

TSY
2-12-09

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 2-12-09 TIME 0652

INTOXILYZER SERIAL # 80-001163 LOCATION DUI HQ DDW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
02/12/2009

06:59:31

2 of 9
On

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

154

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:00:22
Diagnostic Test	Pass	07:00:50
Air Blank	0.000	07:01:17
0.100 Cal Check	0.099	07:01:36
Air Blank	0.000	07:02:05
Subject Test	0.000	07:02:28
Air Blank	ABT*	07:02:34
Air Blank	0.000	07:03:02

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
02/12/2009

07:03:55

3069
Cu
134

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:04:32
Air Blank	IPS**	07:04:42

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
02/12/2009

07:05:27

4 of 9
cu

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:06:16
Diagnostic Test	Pass	07:06:44
Air Blank	0.000	07:07:11
0.100 Cal Check	0.099	07:07:30
Air Blank	0.000	07:07:59
Subject Test	0.441*	07:11:08
Air Blank	0.000	07:11:40
Air Blank	0.000	07:12:07
Subject Test	XXX**	07:12:50
Air Blank	0.000	07:13:23

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
02/12/2009 07:13:48

569
cu
754

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:14:16
Cal Check	0.098	07:14:34
Air Blank	0.000	07:15:04

Clare Conley
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6869
CW
TSY
2-12-09

QA Specialist Clara Conley AGENCY Tucson Police Department
DATE 2-12-09 TIME 0652
Intoxilyzer Serial # 80-001163 Location ODW DW HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clara Conley

DPS Form Exh G-2 (Rev. 05-01)

↓↓↓

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
02/12/2009 06:52:56

7569
C
TS4

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:53:24
Cal Check	0.099	06:53:42
Air Blank	0.000	06:54:12

Clare Conley

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8069
lu

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TSY
2-12-09

Date: 2-12-09
 QAS: Clare Conley
 Intoxilyzer Serial Number: 80-001163
 New Lot Number: 804302I
 New Expiration Date: 2-17-10

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
02/12/2009 07:20:38

9 of 9
Cu

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

154

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:21:05
Cal Check	0.099	07:21:24
Air Blank	0.000	07:21:53

Clare Conley
Operator's Signature