



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
 3-17-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/12/09 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 2/12/09 and 3/12/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —  
 The instrument was returned to service on —

Sarah Cione

Sarah L. Cione


3/17/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17, day of March, 20 09

Ant Hernandez

Notary Public

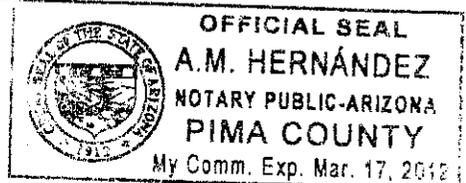

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
3-17-09

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/12/09 TIME 9:34

INTOXILYZER SERIAL # 80-001163 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7sc

154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
03/12/2009

09:37:15

Standard Lot#: 804302I  
Last Changed By: C.. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:38:02
Diagnostic Test	Pass	09:38:29
Air Blank	0.000	09:38:56
0.100 Cal Check	0.101	09:39:15
Air Blank	0.000	09:39:44
Subject Test	0.000	09:40:39
Air Blank	0.000	09:41:06
Five Minute Wait	ABT*	
Air Blank	0.000	09:41:50

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7 sc

154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
03/12/2009 09:42:04

Standard Lot#: 804302I  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:42:36
Air Blank	IPS**	09:42:44

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

4 of 7 sc

TS9

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
03/12/2009 09:43:28

Standard Lot#: 804302I  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:44:31
Diagnostic Test	Pass	09:44:59
Air Blank	0.000	09:45:26
0.100 Cal Check	0.099	09:45:45
Air Blank	0.000	09:46:14
Subject Test	0.168*	09:49:23
Air Blank	0.000	09:49:55
Air Blank	0.000	09:50:23
Subject Test	XXX**	09:50:43
Air Blank	0.000	09:51:17

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7 sc

TS4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
03/12/2009 09:51:41

-----  
Standard Lot#: 804302I  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:52:08
Cal Check	0.099	09:52:27
Air Blank	0.000	09:52:57

Sarah Cione  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSY  
3-17-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/12/09 TIME 9:34  
Intoxilyzer Serial # 80-001163 Location DUT HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
03/12/2009 09:34:00

7 of 7 sc

TS4

-----  
Standard Lot#: 804302I  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:34:27
Cal Check	0.100	09:34:46
Air Blank	0.000	09:35:15



Operator's Signature