



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

*By
2-16-2010*

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 02/11/10 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 01/22/10 and 02/11/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

02/11/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 11th day of February, 2010.

Amy Hernandez

Notary Public



*Includes documents specified in A.A.C. R13-10-104.A.5

187
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TBY
2-16-2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 02/11/10 TIME 11:01

INTOXILYZER SERIAL # 80-001163 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

277
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
02/11/2010

11:04:15

BH

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:05:08
Diagnostic Test	Pass	11:05:36
Air Blank	0.000	11:06:03
0.100 Cal Check	0.096	11:06:22
Air Blank	0.000	11:06:51
Subject Test	0.000	11:07:12
Air Blank	ABT*	11:07:18
Air Blank	0.000	11:07:47

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

317
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
02/11/2010 11:08:05

BA

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:08:47
Air Blank	IPS**	11:08:54

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

487
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
02/11/2010

11:09:11

B4

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:10:02
Diagnostic Test	Pass	11:10:30
Air Blank	0.000	11:10:57
0.100 Cal Check	0.096	11:11:16
Air Blank	0.000	11:11:46
Subject Test	0.043*	11:14:54
Air Blank	0.000	11:15:27
Air Blank	0.000	11:15:54
Subject Test	XXX**	11:16:30
Air Blank	0.000	11:17:04

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

577
CS


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
02/11/2010 11:17:27

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

134

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:17:55
Cal Check	0.095	11:18:13
Air Blank	0.000	11:18:43

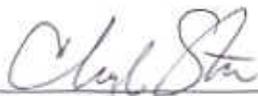

Operator's Signature

EXHIBIT G-2

697
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

By
2-16-2010

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
 DATE 02/11/10 TIME 11:01
 Intoxilyzer Serial # 80-001163 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
 OR
 Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
02/11/2010


11:01:54

Sh

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:02:21
Cal Check	0.096	11:02:40
Air Blank	0.000	11:03:09

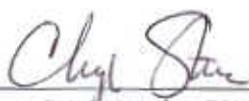

Operator's Signature