



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TSY
3-16-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/16/10 Cal 1 0.096 AC Cal 2 0.097 AC

Between the dates of 2/11/10 and 3/16/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


3/16/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 16th, day of March, 2010.

A.M. Hernandez

Notary Public

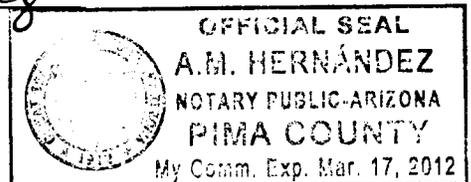

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSY
3.16.2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/16/10 TIME 1507

INTOXILYZER SERIAL # 80-001163 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS 33 days since last standard calibration check procedure

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
03/16/2010

15:09:17

SH

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:10:03
Diagnostic Test	Pass	15:10:30
Air Blank	0.000	15:10:58
0.100 Cal Check	0.097	15:11:16
Air Blank	0.000	15:11:46
Subject Test	0.000	15:12:07
Air Blank	0.000	15:12:35
Five Minute Wait	ABT*	
Air Blank	0.000	15:13:12

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7 sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
03/16/2010 15:13:27

Bh

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:14:13
Diagnostic Test	Pass	15:14:40
Air Blank	0.000	15:15:07
0.100 Cal Check	0.096	15:15:26
Air Blank	0.000	15:15:56
Subject Test	RFI*	15:16:13
Air Blank	IPS**	15:16:22

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc

BH

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
03/16/2010

15:16:42

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:17:25
Diagnostic Test	Pass	15:17:52
Air Blank	0.000	15:18:19
0.100 Cal Check	0.097	15:18:38
Air Blank	0.000	15:19:08
Subject Test	0.087*	15:22:16
Air Blank	0.000	15:22:49
Air Blank	0.000	15:23:16
Subject Test	XXX**	15:23:41
Air Blank	0.000	15:24:15

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
03/16/2010

15:24:28

5 of 7 sc

BY

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:24:55
Cal Check	0.097	15:25:14
Air Blank	0.000	15:25:44

Sarah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

By
3-16-2010

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/16/10 TIME 1507
Intoxilyzer Serial # 80-001163 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
03/16/2010 15:07:50

124

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:08:17
Cal Check	0.096	15:08:36
Air Blank	0.000	15:09:06

Sarah Cione
Operator's Signature