



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSy
4.18.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 04/07/11 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 03/10/2011 and 04/07/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


04/18/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 18th, day of April, 20 11.

A.M. Hernandez

Notary Public

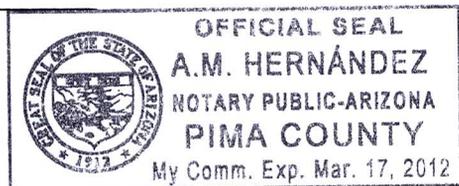

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

197
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

T34
4-18-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl Straub AGENCY Tucson Police Department
DATE 04/07/2011 TIME 06:54
INTOXILYZER SERIAL # 80-001163 LOCATION TEAM 3

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

297
CS

TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
04/07/2011

06:55:35

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:56:26
Diagnostic Test	Pass	06:56:54
Air Blank	0.000	06:57:21
0.100 Cal Check	0.098	06:57:40
Air Blank	0.000	06:58:10
Subject Test	0.000	06:58:30
Air Blank	ABT*	06:58:37
Air Blank	0.000	06:59:05

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

397
CS

TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
04/07/2011

06:59:13

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:59:47
Air Blank	IPS**	06:59:54

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

497
08

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
04/07/2011

06:59:59

TS4

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:00:44
Diagnostic Test	Pass	07:01:11
Air Blank	0.000	07:01:39
0.100 Cal Check	0.098	07:01:58
Air Blank	0.000	07:02:27
Subject Test	0.162*	07:05:36
Air Blank	0.000	07:06:08
Air Blank	0.000	07:06:35
Subject Test	XXX**	07:07:08
Air Blank	0.000	07:07:43

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence


5 of 7
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
04/07/2011 07:08:09

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:08:36
Cal Check	0.097	07:08:55
Air Blank	0.000	07:09:25

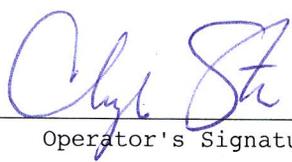
Operator's Signature

697
CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

BG
4.18.2011

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 04/07/2011 TIME 06:54
Intoxilyzer Serial # 80-001163 Location TEAM 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

797
08


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
04/07/2011 06:54:12

134

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:54:40
Cal Check	0.099	06:54:58
Air Blank	0.000	06:55:28

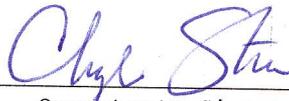

Operator's Signature