



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
1/19/12

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1.12.2012 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 12.16.2011 and 1.12.2012 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on                     .  
The instrument was returned to service on                     .

Terry S. Gallegos

Terry S. Gallegos


1.18.2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 18, day of January, 2012.

A.M. Hernandez

Notary Public

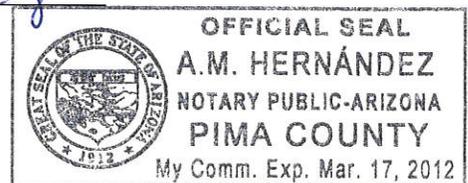

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10/9  
TSG

SC ✓  
1/19/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 1-12-2012 TIME 1051

INTOXILYZER SERIAL # 80-001163 LOCATION ODM

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning. ✓ ✓ ✓  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system. ✓  
Not a Successfully Completed Test Sequence printed ✓  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
01/12/2012

10:53:37

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:54:44
Diagnostic Test	Pass	10:55:11
Air Blank	0.000	10:55:38
0.100 Cal Check	0.096	10:55:57
Air Blank	0.000	10:56:27
Subject Test	0.000	10:57:08
Air Blank	0.000	10:57:35
Five Minute Wait	ABT*	
Air Blank	0.000	10:58:10

\*Sequence Aborted

20 of 9  
TS4  
SC

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
01/12/2012

10:58:15

30f9  
TS4  
SC

---

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:59:09
Diagnostic Test	Pass	10:59:36
Air Blank	RFI*	10:59:53
Air Blank	IPS**	11:00:01

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
01/12/2012

11:00:09

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST,SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:01:03
Diagnostic Test	Pass	11:01:30
Air Blank	0.000	11:01:58
0.100 Cal Check	0.097	11:02:17
Air Blank	0.000	11:02:46
Subject Test	0.035*	11:05:55
Air Blank	0.000	11:06:26
Air Blank	0.000	11:06:54
Subject Test	XXX**	11:07:18
Air Blank	0.000	11:07:53

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

4 of 9  
B4  
SC

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
01/12/2012 11:08:56

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

50F9  
TS4  
SC

Test	g/210L	Time
Air Blank	0.000	11:09:23
Cal Check	0.096	11:09:42
Air Blank	0.000	11:10:12

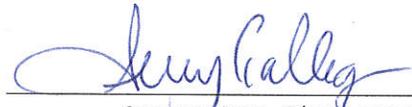
  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 9  
TS4  
SCV  
1/19/12

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 1.12.2012 TIME 1051  
Intoxilyzer Serial # 80-001163 Location ODM

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.


INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
01/12/2012 10:51:02

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

7089  
DL  
SC

Test	g/210L	Time
Air Blank	0.000	10:51:30
Cal Check	0.096	10:51:48
Air Blank	0.000	10:52:18

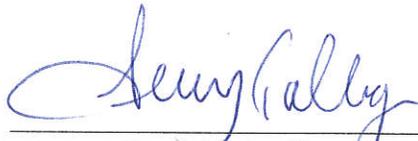  
-----  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

80f.9  
DS4  
SC  
1/19/12

Date: 1.12.2012  
 QAS: Terry Gallegos  
 Intoxilyzer Serial Number: 80-001163  
 New Lot Number: 23511100A2  
 New Expiration Date: 11.1.2013

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
 Select "D" for dry gas and enter the following information:  
     Enter Target Value: 0.100  
     Enter Control Lot #: From the tank label  
     Enter Exp Date: From the tank label (mm/dd/yyyy)  
     Edit info Y/N? Y or N  
     Update Std Info? Y ENTER  
     "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."  
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
01/12/2012 11:13:26


-----  
Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

9 of 9  
TS4

SC

Test	g/210L	Time
Air Blank	0.000	11:13:53
Cal Check	0.099	11:14:12
Air Blank	0.000	11:14:42

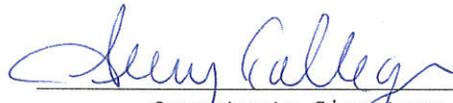  
\_\_\_\_\_  
Operator's Signature