



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 4.11.2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/5/12 Cal 1 0.098 AC Cal 2 0.099 AC

Between the dates of 3/9/12 and 4/5/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/9/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th, day of April, 2012.

Elma Escarcega

Notary Public

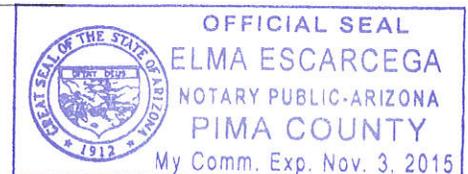

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

T54
4.11.2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/5/12 TIME 1436

INTOXILYZER SERIAL # 80-001163 LOCATION ODM-Tm 3

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7sc

TSG

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
04/05/2012 14:38:21

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:39:09
Diagnostic Test	Pass	14:39:37
Air Blank	0.000	14:40:04
0.100 Cal Check	0.099	14:40:23
Air Blank	0.000	14:40:53
Subject Test	0.000	14:41:17
Air Blank	0.000	14:41:44
Five Minute Wait	ABT*	
Air Blank	0.000	14:42:23

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{sc}

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
04/05/2012

14:42:33

TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:43:19
Diagnostic Test	Pass	14:43:46
Air Blank	RFI*	14:43:55
Air Blank	IPS**	14:44:03

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{sc}

TSG

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
04/05/2012 14:44:14

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:44:57
Diagnostic Test	Pass	14:45:24
Air Blank	0.000	14:45:52
0.100 Cal Check	0.100	14:46:10
Air Blank	0.000	14:46:40
Subject Test	0.106*	14:49:48
Air Blank	0.000	14:50:21
Air Blank	0.000	14:50:48
Subject Test	XXX**	14:51:06
Air Blank	0.000	14:51:40

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
04/05/2012 14:51:54

5 of 7 sc

Bh

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:52:21
Cal Check	0.099	14:52:40
Air Blank	0.000	14:53:10



Operator's Signature

EXHIBIT G-2

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSH
4.11.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 4/5/12 TIME 1436
Intoxilyzer Serial # 80-001163 Location ODM-Tm 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
04/05/2012 14:36:48

7 of 7sc

TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:37:15
Cal Check	0.098	14:37:34
Air Blank	0.000	14:38:04



Operator's Signature