



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

SC  
 2/1/13

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-11-2013 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of — and — the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
 The instrument was returned to service on 2-11-2013.

Terry Gallegos

Terry S. Gallegos


2-26-2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 26<sup>th</sup>, day of February, 2013.

Adma Escarcega

Notary Public

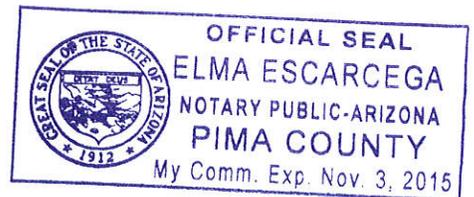

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

10F12  
TJh

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

5c  
3/1/13

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegas 28847 AGENCY Tucson Police Department  
DATE 2-11-2013 TIME 1830  
INTOXILYZER SERIAL # 80-001163 LOCATION Crime Lab

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegas

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
02/11/2013

18:37:17

20f12  
TD4

5c

---

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	18:38:07
Diagnostic Test	Pass	18:38:42
Air Blank	0.000	18:39:10
0.100 Cal Check	0.098	18:39:31
Air Blank	0.000	18:40:00
Subject Test	0.000	18:40:23
Air Blank	0.000	18:40:51
Wait	ABT*	
Air Blank	0.000	18:41:28

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

30f12  
TS4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
02/11/2013 18:41:56

SC

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:42:50
Diagnostic Test	Pass	18:43:25
Air Blank	0.000	18:43:53
0.100 Cal Check	RFI*	18:44:09
Air Blank	IPS**	18:44:17

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
02/11/2013 18:44:24

4 of 12  
TS4

5c

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:45:14
Diagnostic Test	Pass	18:45:49
Air Blank	0.000	18:46:16
0.100 Cal Check	0.097	18:46:37
Air Blank	0.000	18:47:06
Subject Test	0.244*	18:50:14
Air Blank	0.000	18:50:47
Air Blank	0.000	18:51:14
Subject Test	XXX**	18:51:42
Air Blank	0.000	18:52:15

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
02/11/2013 18:53:47

50812  
TB4

-----  
Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847  
-----

30

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	18:54:18
0.100 Cal Check	0.097	18:54:39
Air Blank	0.000	18:55:08

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

60f12  
JSH

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
3/1/13

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 2-11-2013 TIME 1830  
Intoxilyzer Serial # 80-001163 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
02/11/2013 18:30:57

7 of 12  
Dy

-----  
Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847  
-----

5<

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	18:31:24
0.100 Cal Check	0.097	18:31:45
Air Blank	0.000	18:32:15



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

80512

TJH

SC  
3/1/13

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 2.11.2013  
QAS: Terry Gallegos  
Intoxilyzer Serial Number: 80-001163  
New Lot Number: 23511100AZ  
New Expiration Date: 11-1-2013

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove ~~and sign printed~~ record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
02/11/2013 18:29:25

90f12  
B4

-----  
Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847  
-----

SK

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	18:29:53
0.100 Cal Check	0.098	18:30:14
Air Blank	0.000	18:30:43



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

10 of 12  
TB  
SC

INTOXILYZER SERVICE RECORD

QAS: Terry Gallegos 28847

Intoxilyzer Serial Number: 80-001163

Date Returned to Service: 2.11.2013

11 of 12 TS4

5L

# Certificate of Calibration

This is to certify the calibration of **Intoxilyzer**® serial number 80-001163, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to conform to the National Highway Traffic Safety Administration (NHTSA) Standard for Devices to Measure Breath Alcohol (Federal Register, Vol.58, No.179, pp 48705-48710, Sept. 17, 1993) for accuracy and precision. Reference materials are traceable through the National Institute of Standards and Technology (NIST) to the International System of Units (SI).

Date Feb. 4, 2013 Signed Dale Rohrer  
Technician



316 East 9<sup>th</sup> Street  
Owensboro, KY 42303  
USA

Part No. 650517 Rev.A

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

