



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

By
5.13.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/1/13 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 4/3/13 and 5/1/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


5/6/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 6th, day of May, 2013.

Elma Escarcega

Notary Public

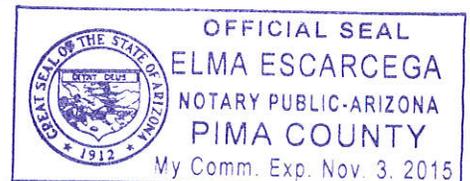

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc
TS4
5.13.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/1/13 TIME 1611

INTOXILYZER SERIAL # 80-001163 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
05/01/2013

16:13:35

2 of 7_{sc}
TJ

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:14:18
Diagnostic Test	Pass	16:14:53
Air Blank	0.000	16:15:21
0.100 Cal Check	0.097	16:15:41
Air Blank	0.000	16:16:11
Subject Test	0.000	16:16:55
Air Blank	0.000	16:17:22
Wait	ABT*	
Air Blank	0.000	16:18:13

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
05/01/2013 16:18:20

3 of 7 sc
TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:19:16
Diagnostic Test	Pass	16:19:51
Air Blank	0.000	16:20:18
0.100 Cal Check	0.097	16:20:39
Air Blank	0.000	16:21:09
Subject Test	RFI*	16:21:27
Air Blank	IPS**	16:21:44

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
05/01/2013 16:21:54

4 of 7 sc
784

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:22:38
Diagnostic Test	Pass	16:23:13
Air Blank	0.000	16:23:41
0.100 Cal Check	0.097	16:24:01
Air Blank	0.000	16:24:31
Subject Test	0.249*	16:27:38
Air Blank	0.000	16:28:10
Air Blank	0.000	16:28:37
Subject Test	XXX**	16:29:01
Air Blank	0.000	16:29:35

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
05/01/2013 16:29:58

5 of 7
sc
TJh

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:30:27
0.100 Cal Check	0.096	16:30:48
Air Blank	0.000	16:31:17

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
sc

TJG

5.13.2013

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 5/1/13 TIME 1611
Intoxilyzer Serial # 80-001163 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
05/01/2013 16:11:35

7 of 7 SC
TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:12:03
0.100 Cal Check	0.096	16:12:24
Air Blank	0.000	16:12:53

Preventative Maintenance
Performed -- Timer Reset