



CERTIFICATION by CUSTODIAN of RECORDS *jk*
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-1-2006 Cal 1 0.100 AC Cal 2 0.101 AC

The instrument was placed into service on 12-1-2006.


Terry S. Gallegos

Terry S. Gallegos

1-4-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.

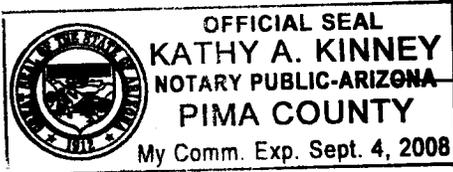

Kathy A. Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

10 of 8
TS4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE *jk*

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 12-1-2006 TIME 1244
INTOXILYZER SERIAL # 80-00408 LOCATION Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning. ✓✓✓
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed ✓
Deficient sample printed. ✓
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE *Terry Gallegos*

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

258
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT *je*
Serial Number : 80-001108
12/01/2006 12:45:44

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:47:00
Diagnostic Test	Pass	12:47:27
Air Blank	0.000	12:47:54
0.100 Cal Check	0.101	12:48:13
Air Blank	0.000	12:48:43
Subject Test	0.000	12:48:59
Air Blank	0.000	12:49:27
Five Minute Wait	ABT*	
Air Blank	0.000	12:50:02

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 8
P4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
12/01/2006 12:50:12

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	12:50:54
Air Blank	0.000	12:51:22

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8
TD4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT *ye*
Serial Number : 80-001108
12/01/2006 12:51:44

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:52:47
Diagnostic Test	Pass	12:53:15
Air Blank	0.000	12:53:42
0.100 Cal Check	0.100	12:54:01
Air Blank	0.000	12:54:31
Subject Test	RFI*	12:54:50
Air Blank	0.000	12:55:18

*RFI Detect

Not a Successfully Completed Test Sequence

5088
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT *je*
Serial Number : 80-001108
12/01/2006 12:55:28

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:56:41
Diagnostic Test	Pass	12:57:09
Air Blank	0.000	12:57:36
0.100 Cal Check	0.100	12:57:55
Air Blank	0.000	12:58:24
Subject Test	0.122*	13:01:33
Air Blank	0.000	13:02:05
Air Blank	0.000	13:02:33
Subject Test	XXX**	13:03:19
Air Blank	0.000	13:03:53

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

6078
TSh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT *je*
Serial Number: 80-001108
12/01/2006 13:05:14

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	13:05:41
Cal Check	0.101	13:06:00
Air Blank	0.000	13:06:30


Operator's Signature

758
754

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)


ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 12.1.2000 TIME 1244
Intoxilyzer Serial # 80-001108 Location Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

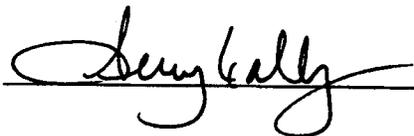SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

8 of 8
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001108
12/01/2006 12:43:23

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	12:43:50
Cal Check	0.100	12:44:09
Air Blank	0.000	12:44:39



Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

Handwritten initials

INTOXILYZER SERVICE RECORD

*1 of 1
R4*

Intoxilyzer Serial Number: 80-001108

Initial Service Date: 12.1.2006

CMI, Inc. Certificate of Calibration Date: 2.4.2005