



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

B4

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/22/06 Cal 1 0.098 AC Cal 2 0.096 AC


Between the dates of 12/1/06 and 12/22/06 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Jennifer L. Coombs
Jennifer L. Coombs

1/4/07
ADPS QAS Permit # 33120

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.

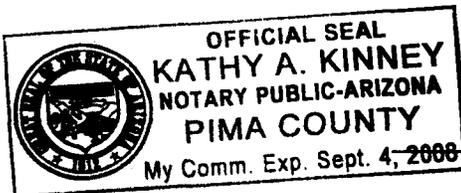

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 9
jcc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES *TS4*
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Jennifer Coombs 46160 AGENCY Tucson Police Department
DATE December 22, 2006 TIME 1545
INTOXILYZER SERIAL # 80-001108 LOCATION DDM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Jennifer L. Coombs

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
12/22/2006

15:56:11

2 of 9
file
TSM

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:57:15
Diagnostic Test	Pass	15:57:43
Air Blank	0.000	15:58:10
0.100 Cal Check	0.098	15:58:29
Air Blank	0.000	15:58:59
Subject Test	0.000	15:59:16
Air Blank	0.000	15:59:43
Five Minute Wait	ABT*	
Air Blank	0.000	16:00:18

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
12/22/2006

16:00:24

BY 3 of 9
ju

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:02:32
Diagnostic Test	Pass	16:02:59
Air Blank	IPS*	16:03:06
Air Blank	0.000	16:03:34

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
12/22/2006

15:50:45

TBh

4 of 9
je

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:52:24
Diagnostic Test	Pass	15:52:51
Air Blank	0.000	15:53:18
0.100 Cal Check	0.098	15:53:37
Air Blank	RFI*	15:53:44
Air Blank	0.000	15:54:12

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
12/22/2006

16:03:52

TS4

5 of 9
ju

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:04:49
Diagnostic Test	Pass	16:05:17
Air Blank	0.000	16:05:44
0.100 Cal Check	0.097	16:06:03
Air Blank	0.000	16:06:33
Subject Test	0.047*	16:09:42
Air Blank	0.000	16:10:14
Air Blank	0.000	16:10:41
Subject Test	XXX**	16:10:58
Air Blank	0.000	16:11:32

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 9
je

ARIZONA DEPARTMENT OF PUBLIC SAFETY

104

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Junifer Coombs 46160 AGENCY Tucson Police Department
DATE December 22, 2006 TIME 1545
Intoxilyzer Serial # 80-001108 Location ODM

✓ ✓
✓ ✓
✓ ✓
✓ ✓
✓ ✓
✓ ✓
✓ ✓

1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.

OR

Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C

2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"

3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.

4. Air blank completed.

5. Calibration check completed. Test results 0. 098 AC. # 0.098AC

6. Air blank completed.

7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Junifer L. Coombs

DPS Form Exh G-2 (Rev. 05-01)

Performed two control tests since paper was placed in the printer upside down for first test.

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

Location: TUCSON POLICE DEPT

Serial Number: 80-001108

12/22/2006

15:44:46

TS

8 of 9
jlc


Standard Lot#: 502602I

Last Changed By: T. GALLEGOS #28847

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	15:45:13
Cal Check	0.098	15:45:32
Air Blank	0.000	15:46:02

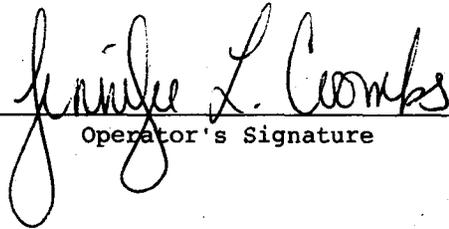
Operator's Signature

