



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

TS4  
7-10-2007

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-6-07 Cal 1 0.097 AC Cal 2 0.098 AC

Between the dates of 6-8-07 and 7-6-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

7-9-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 9, day of July, 2007.

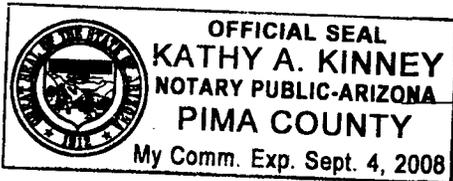

Kathy A Kinney

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 9 m

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

PM

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 7-6-07 TIME 0905

INTOXILYZER SERIAL # 82-001108 LOCATION ODM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2069 cu

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001108  
07/06/2007

09:06:57

PS4

Standard Lot#: 629202I  
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:07:48
Diagnostic Test	Pass	09:08:15
Air Blank	0.000	09:08:42
0.100 Cal Check	0.098	09:09:01
Air Blank	0.000	09:09:31
Subject Test	0.000	09:09:47
Air Blank	0.000	09:10:14
Five Minute Wait	ABT*	
Air Blank	0.000	09:10:48

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

369 cu

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001108  
07/06/2007

09:11:13

TS4

Standard Lot#: 629202I  
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST,SOAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:12:07
Diagnostic Test	Pass	09:12:34
Air Blank	0.000	09:13:02
0.100 Cal Check	RFI*	09:13:15
Air Blank	0.000	09:13:46

\*RFI Detect

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 9 cu

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001108  
07/06/2007

09:13:50

Ph

Standard Lot#: 629202I  
Last Changed By: J. COOMBS #46160  
QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	09:14:26
Air Blank	0.000	09:14:54

\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 9 Cu  
9

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001108  
07/06/2007

09:14:59

Sh

Standard Lot#: 629202I  
Last Changed By: J. COOMBS #46160  
OAS: C. CONLEY #34517  
TUCSON PD  
Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111  
15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:15:55
Diagnostic Test	Pass	09:16:22
Air Blank	0.000	09:16:50
0.100 Cal Check	0.099	09:17:09
Air Blank	0.000	09:17:38
Subject Test	0.144*	09:20:47
Air Blank	0.000	09:21:21
Air Blank	0.000	09:21:48
Subject Test	0.183	09:22:16
Air Blank	0.000	09:22:50
Five Minute Wait	ABT**	
Air Blank	0.000	09:23:23

\*Deficient Sample  
\*\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

6 of 9 cu

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001108  
07/06/2007

09:23:55

PH

Standard Lot#: 629202I  
Last Changed By: J. COOMBS #46160  
QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST,SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:24:51
Diagnostic Test	Pass	09:25:19
Air Blank	0.000	09:25:46
0.100 Cal Check	0.097	09:26:05
Air Blank	0.000	09:26:35
Subject Test	XXX*	09:26:53
Air Blank	0.000	09:27:28

\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7869 Cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001103

07/06/2007

09:29:36

TSh

-----  
Standard Lot#: 629202I

Last Changed By: J. COOMBS #46160

-----  
QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:30:03
Cal Check	0.098	09:30:22
Air Blank	0.000	09:30:52

*Clare Conley*

-----  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

8 of 9 cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

BM

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 7-6-07 TIME 0905  
Intoxilyzer Serial # 80-001108 Location DDM

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

9.89 Cu

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001108  
07/06/2007 09:05:31

-----  
Standard Lot#: 629202I  
Last Changed By: J. COOMBS #46160  
-----

TS4

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:05:58
Cal Check	0.097	09:06:17
Air Blank	0.000	09:06:47

*Clare Conley*

-----  
Operator's Signature