



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

Tsh
6-3-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/20/08 Cal 1 0.103 AC Cal 2 0.102 AC

Between the dates of — and — the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on 5/20/08.


Sarah Cione

Sarah L. Cione

5/20/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 20, day of May, 2008.

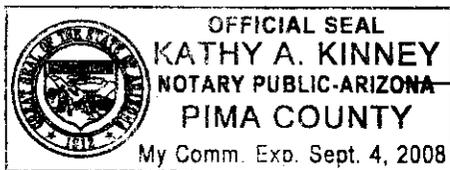

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 12
SC

INTOXILYZER SERVICE RECORD

TS4
6-3-08

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-001108

Date Returned to Service: 5/20/08

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

2 of 12 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
6-3-08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/20/08 TIME 11:47

INTOXILYZER SERIAL # 80-001108 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 102 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 12
SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
05/20/2008

11:49:01

134

Standard Lot#: 725402I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:49:45
Diagnostic Test	Pass	11:50:12
Air Blank	0.000	11:50:40
0.100 Cal Check	0.103	11:50:58
Air Blank	0.000	11:51:28
Subject Test	0.000	11:51:51
Air Blank	0.000	11:52:18
Five Minute Wait	ABT*	
Air Blank	0.000	11:52:55

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
05/20/2008

11:53:41

4 of 12
SC

Standard Lot#: 725402I
Last Changed By: S. CIONE #50616

TSA

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:54:37
Diagnostic Test	Pass	11:55:05
Air Blank	RFI*	11:55:15
Air Blank	IPS**	11:55:23

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

5 of 12
SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
05/20/2008

11:55:52

TSK

Standard Lot#: 725402I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:56:48
Diagnostic Test	Pass	11:57:15
Air Blank	0.000	11:57:43
0.100 Cal Check	0.103	11:58:01
Air Blank	0.000	11:58:31
Subject Test	0.297*	12:01:40
Air Blank	0.000	12:02:13
Air Blank	0.000	12:02:40
Subject Test	XXX**	12:03:06
Air Blank	0.000	12:03:39

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 17, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
05/20/2008

12:04:18

6 of 12
SC

Standard Lot#: 725402I
Last Changed By: S. CIONE #50616

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:04:46
Cal Check	0.102	12:05:05
Air Blank	0.000	12:05:34

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 12
SC

By
6.3.08

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 5/20/08 TIME 11:47
Intoxilyzer Serial # 80-001108 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
05/20/2008

11:47:18

8 of 12
SC
BH

Standard Lot#: 725402I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:47:46
Cal Check	0.103	11:48:04
Air Blank	0.000	11:48:34

Sarah Cione
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

9 of 12
SC

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TBY
6.3.08

Date: 5/20/08
 QAS: Sarah Cione
 Intoxilyzer Serial Number: 80-001108
 New Lot Number: 725402 I
 New Expiration Date: 9/12/2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press
 ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 102 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
05/20/2008

11:45:22

10 of 12
5C

Standard Lot#: 725402I
Last Changed By: S. CIONE #50616

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:45:49
Cal Check	0.102	11:46:08
Air Blank	0.000	11:46:37

Sarah Cione

Operator's Signature

12 of 12
5C
PM

CERTIFICATE OF CALIBRATION

INTOXILYZER® 8000

This is to certify that the calibration of INTOXILYZER serial number 80-001108, manufactured by CMI, subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST (NBS) standard material 1828.

Date 5-08-08 Signed *Dale Howard*
Technician



316 EAST 9TH STREET
OWENSBORO, KY. 42303