



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
8.1.08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/25/08 Cal 1 0.103 AC Cal 2 0.103 AC

Between the dates of 6/27/08 and 7/25/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione


Sarah L. Cione

7/25/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 25th, day of

July, 20 08.

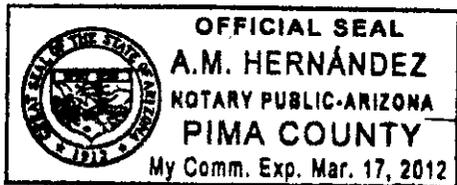

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 9 sc
TSY
8-1-08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/25/08 TIME 9:12

INTOXILYZER SERIAL # 80-001108 LOCATION Team 3

- (✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 103 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 9 sc

184

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
07/25/2008

09:14:10

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:14:53
Diagnostic Test	Pass	09:15:21
Air Blank	0.000	09:15:48
0.100 Cal Check	0.103	09:16:07
Air Blank	0.000	09:16:37
Subject Test	0.000	09:16:58
Air Blank	0.000	09:17:25
Five Minute Wait	ABT*	
Air Blank	0.000	09:18:14

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
07/25/2008

09:18:34

3 of 9 sc
PM

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:19:19
Diagnostic Test	Pass	09:19:46
Air Blank	RFI*	09:19:57
Air Blank	IPS**	09:20:06

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
07/25/2008

09:20:50

4 of 9 SC
184

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:21:36
Diagnostic Test	Pass	09:22:04
Air Blank	0.000	09:22:31
0.100 Cal Check	0.104	09:22:50
Air Blank	0.000	09:23:19
Subject Test	0.078*	09:26:28
Air Blank	0.000	09:26:59
Air Blank	0.000	09:27:27
Subject Test	XXX**	09:27:44
Air Blank	0.000	09:28:16

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
07/25/2008

09:28:46

5 of 9 SC

TS4

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:29:13
Cal Check	0.103	09:29:32
Air Blank	0.000	09:30:01

Sarah Cione

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 9 SC

184
8.1.08

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department

DATE 7/25/08 TIME 9:12

Intoxilyzer Serial # 80-001108 Location Team 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
07/25/2008

09:12:49

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:13:16
Cal Check	0.103	09:13:35
Air Blank	0.000	09:14:05

Sarah Cione

Operator's Signature

7 of 9 SC
184



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 SC

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

104

Date: 7/25/08
 QAS: Sarah Cione
 Intoxilyzer Serial Number: 80-001108
 New Lot Number: 804302 I
 New Expiration Date: 2/17/2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 103 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
07/25/2008 09:11:24

9 of 95C
BM

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:11:52
Cal Check	0.103	09:12:10
Air Blank	0.000	09:12:40

Sarah Cione
Operator's Signature