



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
9-25-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-19-08 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 8-22-08 and 9-19-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


9-25-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 25th, day of September, 2008.

A.M. Hernandez

Notary Public

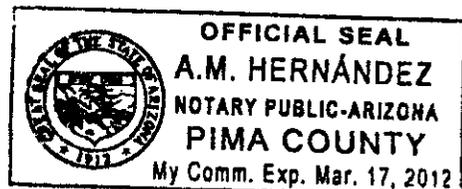

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
in
BY
9.25.08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 9-19-08 TIME 0716

INTOXILYZER SERIAL # 80-001108 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
09/19/2008

07:35:27

2587
lu
Bh

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:36:17
Diagnostic Test	Pass	07:36:44
Air Blank	0.000	07:37:11
0.100 Cal Check	0.102	07:37:30
Air Blank	0.000	07:38:00
Subject Test	0.000	07:38:29
Air Blank	ABT*	07:38:35
Air Blank	0.000	07:39:04

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
09/19/2008

07:39:43

3587
lu
BH

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:40:13
Air Blank	IPS**	07:40:21

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
09/19/2008

07:40:57

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:41:41
Diagnostic Test	Pass	07:42:09
Air Blank	0.000	07:42:36
0.100 Cal Check	0.102	07:42:55
Air Blank	0.000	07:43:24
Subject Test	0.113*	07:46:33
Air Blank	0.000	07:47:04
Air Blank	0.000	07:47:32
Subject Test	XXX**	07:47:51
Air Blank	0.000	07:48:24

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
09/19/2008

07:48:54

5007
TS4 on

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:49:21
Cal Check	0.101	07:49:40
Air Blank	0.000	07:50:09

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

6587
Or
754
9.25.08

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 9-19-08 TIME 0716
Intoxilyzer Serial # 80-001108 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
09/19/2008

07:16:58

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:17:25
Cal Check	0.101	07:17:44
Air Blank	0.000	07:18:13

Clare Conley
Operator's Signature

7567
cu

TS4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 2 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

T54
10-2-2008

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/29/08 TIME 10:23
Intoxilyzer Serial # 80-001108 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 2sc

154

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
09/29/2008 10:23:43

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:24:11
Cal Check	0.102	10:24:30
Air Blank	0.000	10:24:59

Sarah Cione
Operator's Signature