



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

5C  
 11/21/08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11-20-08 Cal 1 0.103 AC Cal 2 0.100 AC

Between the dates of 11-23-08 and 11-20-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
 The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


11-20-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 20<sup>th</sup>, day of November, 2008.

Jim Hernandez

Notary Public

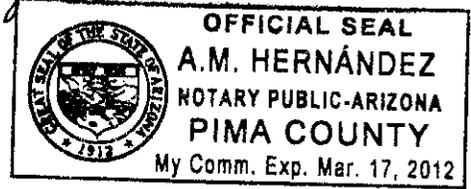

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
11/21/08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 11-20-08 TIME 0850

INTOXILYZER SERIAL # 50-001108 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
11/20/2008

08:52:39

2 of 7  
en

SC

Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:53:32
Diagnostic Test	Pass	08:53:59
Air Blank	0.000	08:54:26
0.100 Cal Check	0.102	08:54:45
Air Blank	0.000	08:55:14
Subject Test	0.000	08:55:36
Air Blank	ABT*	08:55:42
Air Blank	0.000	08:56:11

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
11/20/2008

08:58:10

3-17  
CW  
SC

---

Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	08:58:42
Air Blank	IPS**	08:58:55

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
11/20/2008

08:59:37

4567  
Cu  
SC

Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:00:23
Diagnostic Test	Pass	09:00:51
Air Blank	0.000	09:01:18
0.100 Cal Check	0.102	09:01:37
Air Blank	0.000	09:02:06
Subject Test	0.252*	09:05:14
Air Blank	0.000	09:05:47
Air Blank	0.000	09:06:14
Subject Test	XXX**	09:06:43
Air Blank	0.000	09:07:17

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
11/20/2008

09:07:39

5 of 7  
on  
SC

-----  
Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517  
TUCSON PD -  
-----

Test	g/210L	Time
Air Blank	0.000	09:08:07
Cal Check	0.100	09:08:25
Air Blank	0.000	09:08:55

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6567  
CW

JC  
11/21/08

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 11-20-08 TIME 0850  
Intoxilyzer Serial # 80-001108 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
11/20/2008

08:50:17

-----  
Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:50:44
Cal Check	0.103	08:51:03
Air Blank	0.000	08:51:33

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

7 of 7  
cm

SC