



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY  
3-17-09

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/6/09 Cal 1 0.106 AC Cal 2 0.100 AC

Between the dates of 2/12/09 and 3/6/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —  
The instrument was returned to service on —

Sarah Cione

Sarah L. Cione


3/17/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17, day of March, 20 09.

A.M. Hernandez

Notary Public

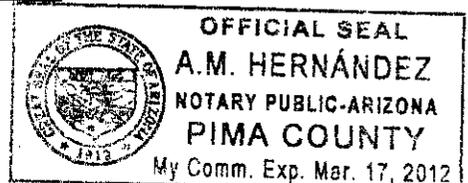

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
3-17-09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 3/6/09 TIME 9:35  
INTOXILYZER SERIAL # 80-001108 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

Replaced the regulator after the first calibration check  
as instrument would not allow an evidentiary test to proceed  
due to gas tank pressure not being recognized.

SIGNATURE

Sarah Cione

2 of 7sc

INTOXILYZER 8000  
 Location: ARIZONA  
 Serial Number : 80-001108  
 03/06/2009 09:55:25

---

Standard Lot#: 804302I  
 Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
 TUCSON PD

---

Operator: S. CIONE #50616  
 TUCSON PD

---

Subject: SQAP, TEST  
 DOB: 11/11/1111  
 Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

156

Test	g/210L	Time
Air Blank	0.000	09:56:26
Diagnostic Test	Pass	09:56:54
Air Blank	0.000	09:57:21
0.100 Cal Check	0.101	09:57:40
Air Blank	0.000	09:58:10
Subject Test	0.000	09:58:29
Air Blank	0.000	09:58:56
Five Minute Wait	ABT*	
Air Blank	0.000	09:59:32

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
03/06/2009

09:59:38

3 of 7sc

154

Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:00:21
Air Blank	IPS**	10:00:30

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
03/06/2009

10:01:18

4 of 7 sc

154

Standard Lot#: 804302I  
Last Changed By: S. CIONE - #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:02:05
Diagnostic Test	Pass	10:02:32
Air Blank	0.000	10:02:59
0.100 Cal Check	0.101	10:03:18
Air Blank	0.000	10:03:48
Subject Test	0.102*	10:06:56
Air Blank	0.000	10:07:28
Air Blank	0.000	10:07:55
Subject Test	XXX**	10:08:28
Air Blank	0.000	10:09:00

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
03/06/2009

10:09:45

5 of 7sc

-----  
Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616  
-----

124

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:10:13
Cal Check	0.100	10:10:32
Air Blank	0.000	10:11:01

*Sarah Cione*

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

TSH  
3-17-09

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/6/09 TIME 9:35  
Intoxilyzer Serial # 80-001108 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 106 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
03/06/2009

09:35:34

7 of 7 sc

-----  
Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616  
-----

154

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:36:01
Cal Check	0.106	09:36:20
Air Blank	0.000	09:36:50

*Sarah Cione*

Operator's Signature