



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

TSY  
4.13.09

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/10/09 Cal 1 0.103 AC Cal 2 0.102 AC

Between the dates of N/A and N/A the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on 4/10/09.

Sarah Cione

Sarah L. Cione


4/13/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 13<sup>th</sup>, day of April, 2009.

Art Hernandez

Notary Public

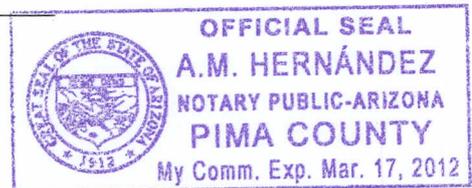

\*Includes documents specified in A.A.C. R13-10-104.A.5



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 10<sub>SC</sub>

INTOXILYZER SERVICE RECORD

TS4  
4.13.09

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-001108

Date Returned to Service: 4/10/09

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

2 of 10 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
4.13.09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/10/09 TIME 15:10

INTOXILYZER SERIAL # 80-001108 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 102 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
04/10/2009

15:13:41

3 of 10<sub>sc</sub>  
124

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:14:29
Diagnostic Test	Pass	15:14:56
Air Blank	0.000	15:15:23
0.100 Cal Check	0.104	15:15:42
Air Blank	0.000	15:16:12
Subject Test	0.000	15:16:34
Air Blank	0.000	15:17:01
Five Minute Wait	ABT*	
Air Blank	0.000	15:17:36

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
04/10/2009

15:17:42

4 of 10<sub>sc</sub>

BM

---

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	15:18:11
Air Blank	IPS**	15:18:19

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

5 of 10sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
04/10/2009

15:18:24

PM

Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:19:10
Diagnostic Test	Pass	15:19:37
Air Blank	0.000	15:20:04
0.100 Cal Check	0.104	15:20:23
Air Blank	0.000	15:20:53
Subject Test	0.101*	15:24:02
Air Blank	0.000	15:24:34
Air Blank	0.000	15:25:01
Subject Test	XXX**	15:25:25
Air Blank	0.000	15:25:58

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
04/10/2009

15:26:39

6 of 10 sc

BY

-----  
Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:27:06
Cal Check	0.102	15:27:25
Air Blank	0.000	15:27:55

*Sarah Cione*

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

7 of 10<sub>sc</sub>

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

Tsh  
4.13.09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 4/10/09 TIME 15:10  
Intoxilyzer Serial # 80-001108 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
04/10/2009

15:10:58

8 of 10 sc

-----  
Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616  
-----

TSY

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:11:26
Cal Check	0.103	15:11:44
Air Blank	0.000	15:12:14

Sarah Cione  
Operator's Signature



TUCSON POLICE DEPARTMENT 9 of 10 sc  
 CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
 FOR  
 CALIBRATION GAS TANK CHANGE

Date: 4/10/09 1544.13.09  
 QAS: Sarah Cione  
 Intoxilyzer Serial Number: 80-001108  
 New Lot Number: 826001 I  
 New Expiration Date: 9/24/2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
 Select "D" for dry gas and enter the following information:  
     Enter Target Value: 0.100  
     Enter Control Lot #: From the tank label  
     Enter Exp Date: From the tank label (mm/dd/yyyy)  
     Edit info Y/N? Y or N  
     Update Std Info? Y ENTER  
     "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 105 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
04/10/2009

15:08:35

-----  
Standard Lot#: 826001I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:09:02
Cal Check	0.105	15:09:21
Air Blank	0.000	15:09:51

*Sarah Cione*  
-----  
Operator's Signature

10 of 10<sub>sc</sub>

TS4