



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

By
8.23.2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/19/10 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 7/22/10 and 8/19/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


8/23/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 23, day of August, 20 10.

A.M. Hernandez

Notary Public

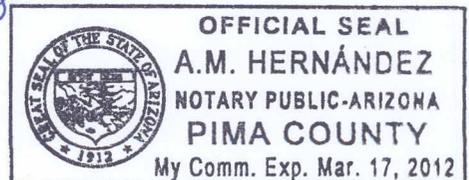

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

DY
8.23.2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 8/19/10 TIME 857
INTOXILYZER SERIAL # 80-001108 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7 sc

PA

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
08/19/2010 08:58:47

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:59:35
Diagnostic Test	Pass	09:00:02
Air Blank	0.000	09:00:30
0.100 Cal Check	0.101	09:00:48
Air Blank	0.000	09:01:18
Subject Test	0.000	09:01:59
Air Blank	0.000	09:02:27
Five Minute Wait	ABT*	
Air Blank	0.000	09:03:04

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
08/19/2010 09:03:38

3 of 7 sc
124

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:04:08
Air Blank	IPS**	09:04:16

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc

SPH

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
08/19/2010 09:04:55

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:05:39
Diagnostic Test	Pass	09:06:07
Air Blank	0.000	09:06:34
0.100 Cal Check	0.101	09:06:53
Air Blank	0.000	09:07:22
Subject Test	0.130*	09:10:31
Air Blank	0.000	09:11:04
Air Blank	0.000	09:11:32
Subject Test	XXX**	09:12:13
Air Blank	0.000	09:12:46

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
08/19/2010

09:13:25

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:13:53
Cal Check	0.100	09:14:12
Air Blank	0.000	09:14:41

Sarah Cione

Operator's Signature

5 of 7 sc

104

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

By
8-23-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 8/19/10 TIME 857
Intoxilyzer Serial # 80-001108 Location DUT HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
08/19/2010

08:57:05

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:57:32
Cal Check	0.101	08:57:51
Air Blank	0.000	08:58:21

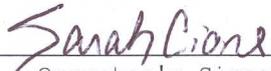

Operator's Signature

7 of 7sc

SM