



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

SC
2/8/08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-8-08 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 1-10-08 and 2-8-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A
The instrument was returned to service on N/A

Clare Conley

Clare L. Conley


2-8-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 8 day of February, 2008.

Kathy A Kinney

Notary Public

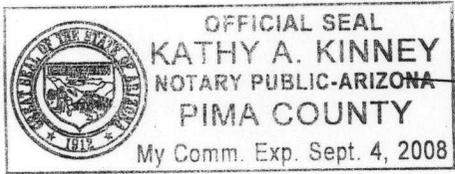

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8 cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SC
2/8/08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 2-8-08 TIME 07:45

INTOXILYZER SERIAL # 80-001165 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
02/08/2008

07:48:47

2/8/08

SC
2/8/08

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:49:37
Diagnostic Test	Pass	07:50:05
Air Blank	0.000	07:50:32
0.100 Cal Check	0.099	07:50:51
Air Blank	0.000	07:51:21
Subject Test	0.000	07:51:43
Air Blank	ABT*	07:51:49
Air Blank	0.000	07:52:18

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
02/08/2008

07:52:26

3/8
cu

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

SC
2/8/08

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:53:00
Air Blank	IPS**	07:53:08

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
02/08/2008

07:53:27

4 of 8
cu

SC
2/8/08

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:54:17
Diagnostic Test	Pass	07:54:45
Air Blank	0.000	07:55:12
0.100 Cal Check	0.099	07:55:31
Air Blank	0.000	07:56:01
Subject Test	XXX*	07:58:35
Air Blank	0.000	07:59:11

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
02/08/2008 08:00:02

5/8
on
SC
2/8/08

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:00:56
Diagnostic Test	Pass	08:01:24
Air Blank	0.000	08:01:51
0.100 Cal Check	0.098	08:02:10
Air Blank	0.000	08:02:41
Subject Test	0.315*	08:05:49
Air Blank	ABT**	08:05:58
Air Blank	0.000	08:06:26

*Deficient Sample
**Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

02/08/2008

08:06:43

Standard Lot#: 725402I

Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:07:10
Cal Check	0.097	08:07:29
Air Blank	0.000	08:07:59

Clare Conley

Operator's Signature

bfg
cu

sc
2/8/08

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7/8
CW

SC
2/8/08

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 2-8-08 TIME 09:45
Intoxilyzer Serial # SD-001165 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
02/08/2008 07:45:55

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

8/8
C
SC
2/8/08

Test	g/210L	Time
Air Blank	0.000	07:46:23
Cal Check	0.099	07:46:42
Air Blank	0.000	07:47:11

C. Conley

Operator's Signature