



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY
 3-15-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 03/10/2011 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 02/10/2011 and 03/10/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
 The instrument was returned to service on N/A.

Cheryl K. Straub
 Cheryl K. Straub


03/15/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 15, day of March, 2011.

A.M. Hernandez

Notary Public

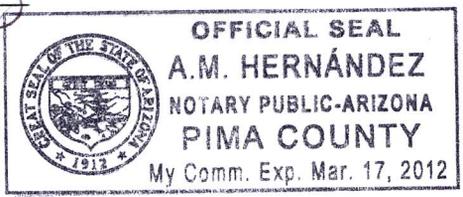

*Includes documents specified in A.A.C. R13-10-104.A.5

187
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
3-15-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 03/10/2011 TIME 07:03

INTOXILYZER SERIAL # 80-001165 LOCATION TEAM 1

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

277
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
03/10/2011 07:04:41

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:05:32
Diagnostic Test	Pass	07:06:00
Air Blank	0.000	07:06:27
0.100 Cal Check	0.098	07:06:46
Air Blank	0.000	07:07:16
Subject Test	0.000	07:07:35
Air Blank	ABT*	07:07:41
Air Blank	0.000	07:08:09

*Sequence Aborted

~~Not a Successfully Completed Test Sequence~~

377
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
03/10/2011 07:08:13

TS4

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:08:45
Air Blank	IPS**	07:08:56

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

437
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
03/10/2011 07:09:00

TSK

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:09:48
Diagnostic Test	Pass	07:10:16
Air Blank	0.000	07:10:43
0.100 Cal Check	0.098	07:11:02
Air Blank	0.000	07:11:32
Subject Test	0.431*	07:14:41
Air Blank	0.000	07:15:15
Air Blank	0.000	07:15:42
Subject Test	XXX**	07:16:12
Air Blank	0.000	07:16:46

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

507
CG

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
03/10/2011 07:17:01

TS4

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:17:29
Cal Check	0.098	07:17:47
Air Blank	0.000	07:18:18



Operator's Signature

EXHIBIT G-2

687
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

By
3.15.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl Straub AGENCY Tucson Police Department
DATE 03/10/2011 TIME 07:03
Intoxilyzer Serial # 80-001165 Location TEAM 1

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

797
CS

124

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
03/10/2011 07:03:20

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:03:47
Cal Check	0.098	07:04:06
Air Blank	0.000	07:04:36



Operator's Signature