



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
5.10.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 05/05/11 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 04/07/2011 and 05/05/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.

The instrument was returned to service on N/A.

Cheryl K. Straub


Cheryl K. Straub

05/10/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 10, day of

May, 20 11.

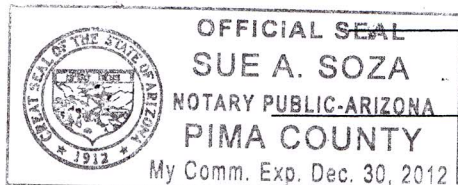

Sue A. Soza

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

187
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

B4
5.10.2011QA SPECIALIST Cheryl Straub AGENCY Tucson Police DepartmentDATE 05/05/2011 TIME 07:18INTOXILYZER SERIAL # 80-001165 LOCATION TEAM I

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
-
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.
- 000
- AC.
-
- (✓) 2. Error recognition logic system functioning.
-
- Not a Successfully Completed Test Sequence printed
-
- (✓) 3. Proper sample recognition system.
-
- Not a Successfully Completed Test Sequence printed
-
- Deficient sample printed.
-
- (✓) 4. Standard Calibration Check standard 0.
- 100
- AC
-
- Result: 0.
- 097
- AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Cheryl Straub

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

05/05/2011

07:20:10

Standard Lot#: 719919

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:20:57
Diagnostic Test	Pass	07:21:25
Air Blank	0.000	07:21:52
0.100 Cal Check	0.098	07:22:11
Air Blank	0.000	07:22:41
Subject Test	0.000	07:22:56
Air Blank	ABT*	07:23:03
Air Blank	0.000	07:23:31

*Sequence Aborted

Not a Successfully Completed Test Sequence

387
CS

BSH

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

05/05/2011

07:23:37

Standard Lot#: 719919

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:24:09
Air Blank	IPS**	07:24:23

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

497
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
05/05/2011

07:24:27

134

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:25:17
Diagnostic Test	Pass	07:25:45
Air Blank	0.000	07:26:12
0.100 Cal Check	0.098	07:26:31
Air Blank	0.000	07:27:01
Subject Test	0.136*	07:30:10
Air Blank	0.000	07:30:42
Air Blank	0.000	07:31:09
Subject Test	XXX**	07:31:31
Air Blank	0.000	07:32:05

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

577
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

05/05/2011

07:32:34

154

Standard Lot#: 719919

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:33:01
Cal Check	0.097	07:33:20
Air Blank	0.000	07:33:50



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl Straub AGENCY Tucson Police Department
DATE 05/05/2011 TIME 07:18
Intoxilyzer Serial # 80-001165 Location TEAM I

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 0.98 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

727
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

05/05/2011

07:18:44

B4


Standard Lot#: 719919

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:19:12
Cal Check	0.098	07:19:31
Air Blank	0.000	07:20:00

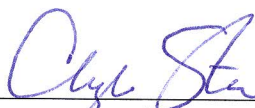

Operator's Signature