



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

SC
11/21/11

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/17/11 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 10/21/11 and 11/17/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


11/21/11

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 21, day of November, 20 11.

A.M. Hernandez

Notary Public

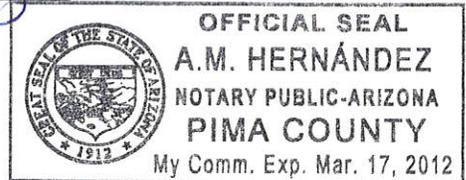

*Includes documents specified in A.A.C. R13-10-104.A.5

18709

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

SC
11/21/11

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST CHERYL STRAUB AGENCY Tucson Police Department

DATE 11/17/11 TIME 07:21

INTOXILYZER SERIAL # 80-001165 LOCATION TEAM 1

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

287cs

SC

INTOXILYZER 8000
 Location: TUCSON POLICE DEPT
 Serial Number : 80-001165
 11/17/2011 07:23:06

Standard Lot#: 04911100A1
 Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
 TUCSON PD

Operator: C. STRAUB #50772
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1111
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:23:52
Diagnostic Test	Pass	07:24:19
Air Blank	0.000	07:24:46
0.100 Cal Check	0.096	07:25:05
Air Blank	0.000	07:25:36
Subject Test	0.000	07:25:52
Air Blank	ABT*	07:25:59
Air Blank	0.000	07:26:28

*Sequence Aborted

Not a Successfully Completed Test Sequence

387
C8

SC

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
11/17/2011 07:26:32

Standard Lot#: 04911100A1
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:27:02
Air Blank	IPS**	07:27:09

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

427
c8

SC

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
11/17/2011 07:27:14

Standard Lot#: 04911100A1
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:27:59
Diagnostic Test	Pass	07:28:27
Air Blank	0.000	07:28:54
0.100 Cal Check	0.096	07:29:13
Air Blank	0.000	07:29:43
Subject Test	0.223*	07:32:52
Air Blank	0.000	07:33:25
Air Blank	0.000	07:33:52
Subject Test	XXX**	07:34:19
Air Blank	0.000	07:34:53

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence


587
08
5C

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
11/17/2011 07:35:11

Standard Lot#: 04911100A1
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:35:39
Cal Check	0.095	07:35:57
Air Blank	0.000	07:36:28

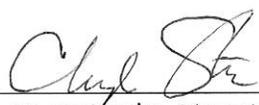

Operator's Signature

EXHIBIT G-2

677
08
SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 11/17/11 TIME 07:21
Intoxilyzer Serial # 80-001165 Location TEAM 1

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

727
08


SC

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
11/17/2011 07:21:39

Standard Lot#: 04911100A1
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:22:06
Cal Check	0.096	07:22:25
Air Blank	0.000	07:22:55

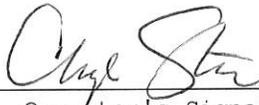

Operator's Signature