



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
12-21-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number

80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-16-11 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 11-17-11 and 12-16-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.

The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

12-16-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 16, day of December, 20 11.

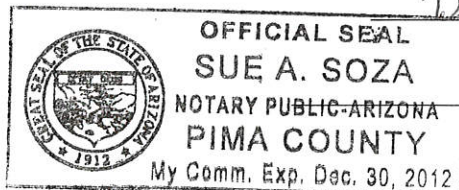

Sue A. Soza

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 12-16-11 TIME 0643

INTOXILYZER SERIAL # 80-001165 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS Added five minutes to time on internal clock.

SIGNATURE Clare Conley

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

12/16/2011

06:45:00

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB

#50772

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:45:44
Diagnostic Test	Pass	06:46:12
Air Blank	0.000	06:46:39
0.100 Cal Check	0.097	06:46:58
Air Blank	0.000	06:47:28
Subject Test	0.000	06:47:49
Air Blank	ABT*	06:47:56
Air Blank	RFI**	06:48:03

*Sequence Aborted

**RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

12/16/2011

06:48:08

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB

#50772

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	06:48:33
Air Blank	0.000	06:49:01

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

12/16/2011

06:49:05

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB

#50772

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:49:52
Diagnostic Test	Pass	06:50:19
Air Blank	0.000	06:50:47
0.100 Cal Check	0.096	06:51:05
Air Blank	0.000	06:51:36
Subject Test	0.440*	06:54:44
Air Blank	0.000	06:55:17
Air Blank	0.000	06:55:44
Subject Test	XXX**	06:56:08
Air Blank	0.000	06:56:41

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

12/16/2011

06:57:31

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:57:58
Cal Check	0.096	06:58:17
Air Blank	0.000	06:58:47



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 12-16-11 TIME 0643
Intoxilyzer Serial # 80-001165 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

12/16/2011

06:43:39

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:44:06
Cal Check	0.097	06:44:25
Air Blank	0.000	06:44:55



Operator's Signature