



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
6-12-2012

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/10/12 Cal 1 0.097 AC Cal 2 0.095 AC

Between the dates of 5/1/12 and 5/10/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 5/10/12.
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


5/24/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 24th day of May, 2012.

Elma Escarcega

Notary Public

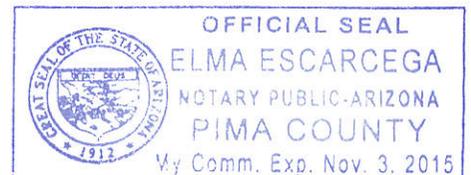

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 9 sc

By
6-12-2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/10/12 TIME 1252

INTOXILYZER SERIAL # 80-001165 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 9 sc
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
05/10/2012 12:54:06

Standard Lot#: 04911100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:55:08
Diagnostic Test	Pass	12:55:36
Air Blank	0.000	12:56:03
0.100 Cal Check	0.097	12:56:22
Air Blank	0.000	12:56:52
Subject Test	0.000	12:57:08
Air Blank	0.000	12:57:36
Five Minute Wait	ABT*	
Air Blank	0.000	12:58:13

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
05/10/2012 12:58:17

3 of 9 sc
Bh

Standard Lot#: 04911100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:59:17
Diagnostic Test	Pass	12:59:44
Air Blank	RFI*	12:59:55
Air Blank	IPS**	13:00:03

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 9_{sc}
BH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
05/10/2012 13:00:12

Standard Lot#: 04911100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:00:55
Diagnostic Test	Pass	13:01:22
Air Blank	0.000	13:01:49
0.100 Cal Check	0.095	13:02:08
Air Blank	0.000	13:02:39
Subject Test	0.102*	13:05:47
Air Blank	0.000	13:06:19
Air Blank	0.000	13:06:46
Subject Test	XXX**	13:07:04
Air Blank	0.000	13:07:37

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence


INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
05/10/2012 13:08:06

5 of 9 sc
Bh

Standard Lot#: 04911100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:08:33
Cal Check	0.095	13:08:52
Air Blank	0.000	13:09:23

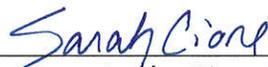

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 9
SC
Bh
6-12-2012

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 5/10/12 TIME 1252
Intoxilyzer Serial # 80-001165 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
05/10/2012 12:52:45

7 of 9 sc
DL

Standard Lot#: 04911100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:53:13
Cal Check	0.097	12:53:31
Air Blank	0.000	12:54:02

Sarah Cione

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 sc
154

INTOXILYZER SERVICE RECORD

Date: 5/10/12

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-001165

Reason For Service: Need to apply more pressure to achieve tone

Date Out of Service: 5/10/12

Calibration Check completed prior to removal from service

Date: 5/10/12

SQAP completed prior to removal from service

Date: 5/10/12

Corrective Action:

Sent to CMI, Inc.

9 of 9 sc

104



Service Evaluation Form

This form **MUST** be completed and enclosed with instrument to be serviced.
Failure to complete and return this form may cause delays in service.

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)
 Name Terry Gallegos Phone: (520) 837-7677
 Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:
Finance / Accounts Payable Tucson Police Department Crime Lab
Po Box 27450 1306 W. Miracle Mile
Tucson, AZ 85726-7450 Tucson, AZ 85705

3. Serial Number: 80-001165 Instrument Model: Intox P000

4. Detailed Description of Problem:
need more pressure to get tone, Static Reading ~70

Hazardous Material Warning! - DO NOT return gas cylinders with instrument!

I Authorize Repairs Up To: All \$250 \$500 \$750 Other \$ _____
 Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:
Terry Gallegos Toxicology Coordinator
 Name (Please Print) Title
Terry Gallegos 5.10.2012
 Signature Date

Ship to:
 CMI, Inc.
 Attn: Service Dept.
 316 East Ninth Street
 Owensboro, KY 42303

No, please send estimate before repairs are made.
 Note: An estimate will be faxed before performing any repairs and may cause delays in service.
 An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.