



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
9-20-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/13/12 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 8/16/12 and 9/13/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

9/17/12


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17th, day of September, 2012.

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

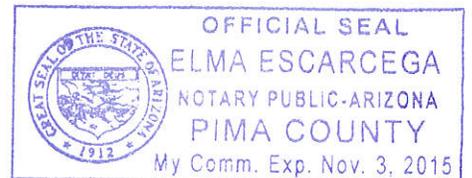

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 7_{sc}

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TS4
9.20.2012

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/13/12 TIME 1205

INTOXILYZER SERIAL # 80-001165 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7sc

TS4

INTOXILYZER 8000
 Location: TUCSON POLICE DEPT
 Serial Number: 80-001165
 Core Version: 8105.48
 09/13/2012 12:07:36

Standard Lot#: 23511100A2
 Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1980
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:08:23
Diagnostic Test	Pass	12:08:58
Air Blank	0.000	12:09:25
0.100 Cal Check	0.102	12:09:46
Air Blank	0.000	12:10:16
Subject Test	0.000	12:10:49
Air Blank	0.000	12:11:17
Wait	ABT*	
Air Blank	0.000	12:11:58

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
09/13/2012 12:12:19

3 of 7sc

TJ4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:13:09
Diagnostic Test	Pass	12:13:45
Air Blank	0.000	12:14:12
0.100 Cal Check	0.101	12:14:33
Air Blank	0.000	12:15:02
Subject Test	RFI*	12:15:29
Air Blank	IPS**	12:15:36

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{sc}

TSh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
09/13/2012 12:15:40

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:16:21
Diagnostic Test	Pass	12:16:56
Air Blank	0.000	12:17:23
0.100 Cal Check	0.101	12:17:44
Air Blank	0.000	12:18:13
Subject Test	0.249*	12:21:21
Air Blank	0.000	12:21:58
Air Blank	0.000	12:22:26
Subject Test	XXX**	12:22:39
Air Blank	0.000	12:23:13

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7_{sc}

D4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
09/13/2012 12:23:28

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:23:57
0.100 Cal Check	0.099	12:24:18
Air Blank	0.000	12:24:47

Preventative Maintenance
Performed -- Timer Reset

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

Bh
9.20.2012

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/13/12 TIME 1205
Intoxilyzer Serial # 80-001165 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- () 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
09/13/2012 12:05:40

7 of 7 sc
TSh

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:06:09
0.100 Cal Check	0.101	12:06:30
Air Blank	0.000	12:06:59

Preventative Maintenance
Performed -- Timer Reset