



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 1.7.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/6/12 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 11/8/12 and 12/6/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

12/17/12


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17th, day of December, 2012.

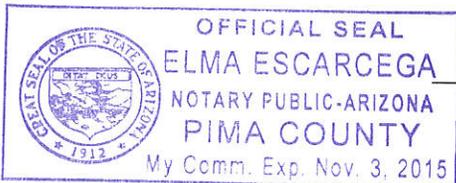

Elma Escarcega

Notary Public

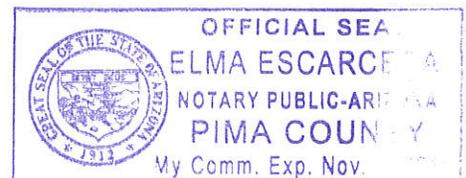

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7sc
B4
1.7.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 12/6/12 TIME 1024

INTOXILYZER SERIAL # 80-001165 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7 sc

134

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/06/2012 10:26:44

Standard Lot#: 13712100A2
Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:27:35
Diagnostic Test	Pass	10:28:09
Air Blank	0.000	10:28:37
0.100 Cal Check	0.101	10:28:58
Air Blank	0.000	10:29:27
Subject Test	0.000	10:29:48
Air Blank	0.000	10:30:16
Wait	ABT*	
Air Blank	0.000	10:30:51

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7 sc
TSK

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/06/2012 10:31:07

Standard Lot#: 13712100A2
Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:31:51
Diagnostic Test	Pass	10:32:26
Air Blank	0.000	10:32:53
0.100 Cal Check	0.101	10:33:14
Air Blank	RFI*	10:33:28
Air Blank	IPS**	10:33:37

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{sc}
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/06/2012 10:38:30

Standard Lot#: 13712100A2
Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:39:19
Diagnostic Test	Pass	10:39:54
Air Blank	0.000	10:40:22
0.100 Cal Check	0.102	10:40:42
Air Blank	0.000	10:41:12
Subject Test	0.262*	10:44:20
Air Blank	0.000	10:44:53
Air Blank	0.000	10:45:20
Subject Test	XXX**	10:45:36
Air Blank	0.000	10:46:10

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7 sc
TB

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/06/2012 10:46:27

Standard Lot#: 13712100A2
Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:46:55
0.100 Cal Check	0.100	10:47:16
Air Blank	0.000	10:47:46

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc
TS4
1-7-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 12/6/12 TIME 1024
Intoxilyzer Serial # 80-001165 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/06/2012 10:24:27

7 of 7 sc
TS4

Standard Lot#: 13712100A2
Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:24:57
0.100 Cal Check	0.101	10:25:18
Air Blank	0.000	10:25:47

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2 SC
B4
1.7.2013

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 11-17-12
 QAS: Off. Honomichuk #42992
 Intoxilyzer Serial Number: 80-001165
 New Lot Number: #13712100A2
 New Expiration Date: 07-01-2014

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 102 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE [Signature]

2 of 2 sc
154

INTOXILYZER 8000
Instrument Initialization
18:31:08 11/17/2012

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
11/17/2012 19:04:23

Standard Lot#: 13712100A2
Last Changed By: W. HONOMICHL #42992

OAS: S. CIONE #50616
TUCSON PD

Operator: W. HONOMICHL #42992
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	19:04:54
0.100 Cal Check	0.102	19:05:14
Air Blank	0.000	19:05:44

Preventative Maintenance NOT
Performed -- Timer NOT Reset