



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TSG  
3.13.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/28/13 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 1/31/13 and 2/28/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

3/5/13


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 5<sup>th</sup>, day of March, 2013.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

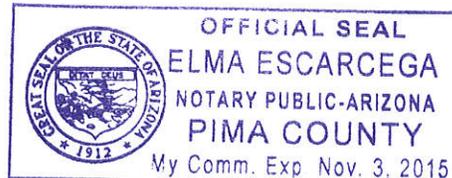

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TS4  
3.13.2013

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/28/13 TIME 1204

INTOXILYZER SERIAL # 80-001165 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7sc  
TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
02/28/2013 12:06:11

---

Standard Lot#: 13712100A2  
Last Changed By: W. HONOMICHL #42992

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:06:57
Diagnostic Test	Pass	12:07:32
Air Blank	0.000	12:08:00
0.100 Cal Check	0.102	12:08:20
Air Blank	0.000	12:08:50
Subject Test	0.000	12:09:08
Air Blank	0.000	12:09:36
Wait	ABT*	
Air Blank	0.000	12:10:12

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7 sc

TS4

INTOXILYZER 8000  
 Location: TUCSON POLICE DEPT  
 Serial Number: 80-001165  
 Core Version: 8105.48  
 02/28/2013 12:15:16

Standard Lot#: 13712100A2  
 Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616  
 TUCSON PD

Operator: S. CIONE #50616  
 TUCSON PD

Subject: SQAP, SQAP, SQAP  
 DOB: 11/11/1980  
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:16:20
Diagnostic Test	Pass	12:16:55
Air Blank	0.000	12:17:22
0.100 Cal Check	0.101	12:17:43
Air Blank	0.000	12:18:12
Subject Test	RFI*	12:18:51
Air Blank	IPS**	12:18:59

\*RFI Detect  
 \*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
02/28/2013

12:20:49

4 of 7 SC  
TS4

---

Standard Lot#: 13712100A2  
Last Changed By: W. HONOMICHL #42992

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:21:47
Diagnostic Test	Pass	12:22:22
Air Blank	0.000	12:22:49
0.100 Cal Check	0.102	12:23:10
Air Blank	0.000	12:23:40
Subject Test	0.169*	12:26:47
Air Blank	0.000	12:27:19
Air Blank	0.000	12:27:47
Subject Test	XXX**	12:28:22
Air Blank	0.000	12:28:57

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
02/28/2013 12:30:53

5 of 7 sc  
TS4

-----  
Standard Lot#: 13712100A2  
Last Changed By: W. HONOMICHL #42992  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:31:21
0.100 Cal Check	0.101	12:31:42
Air Blank	0.000	12:32:11

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 750

TSA

3-13-2013 TSA

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 2/28/13 TIME 1204  
Intoxilyzer Serial # 80-001165 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 102 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
02/28/2013 12:04:41

7 of 7 sc  
BH

-----  
Standard Lot#: 13712100A2  
Last Changed By: W. HONOMICHL #42992  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:05:09
0.100 Cal Check	0.102	12:05:30
Air Blank	0.000	12:06:00

Preventative Maintenance  
Performed -- Timer Reset