



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

TS4  
12-13-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/19/07 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 10/19/07 and 11/19/07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

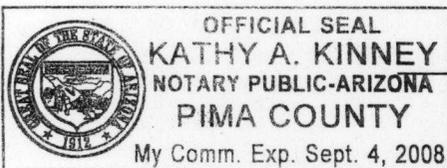
Sarah L. Cione

Sarah L. Cione

11/29/07

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 29, day of November, 2007.



Kathy A Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
12.13.07

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/19/2007 TIME 10:41 am

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS none  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
11/19/2007

10:43:12

2 of 75C

154

Standard Lot#: 719801I  
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:44:04
Diagnostic Test	Pass	10:44:32
Air Blank	0.000	10:44:59
0.100 Cal Check	0.098	10:45:18
Air Blank	0.000	10:45:47
Subject Test	0.000	10:46:07
Air Blank	0.000	10:46:35
Five Minute Wait	ABT*	
Air Blank	0.000	10:47:14

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
11/19/2007

10:47:21

3 of 7sc

154

---

Standard Lot#: 719801I  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:48:11
Diagnostic Test	Pass	10:48:38
Air Blank	RFI*	10:48:47
Air Blank	IPS**	10:49:00

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 26-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
11/19/2007

10:51:11

4 of 7 sc

154

---

Standard Lot#: 719801I  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:52:02
Diagnostic Test	Pass	10:52:30
Air Blank	0.000	10:52:57
0.100 Cal Check	0.098	10:53:16
Air Blank	0.000	10:53:45
Subject Test	0.310*	10:56:54
Air Blank	0.000	10:57:25
Air Blank	0.000	10:57:53
Subject Test	XXX**	10:58:18
Air Blank	0.000	10:58:51

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
11/19/2007 10:59:29

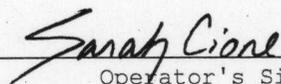
5 of 7 sc

-----  
Standard Lot#: 719801I  
Last Changed By: T. GALLEGOS #28847  
-----

104

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:59:56
Cal Check	0.097	11:00:15
Air Blank	0.000	11:00:44

  
-----  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 sc

TS4

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 11/19/2007 TIME 10:41am  
Intoxilyzer Serial # 80-001199 Location DUT HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
11/19/2007 10:41:29

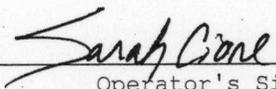
7 of 7 sc

-----  
Standard Lot#: 719801I  
Last Changed By: T. GALLEGOS #28847  
-----

TS4

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:41:57
Cal Check	0.097	10:42:15
Air Blank	0.000	10:42:45

  
-----  
Operator's Signature