



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

Ch Ch
2-20-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/14/08 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 1/16/08 and 2/14/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

2/15/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 15, day of February, 2008.

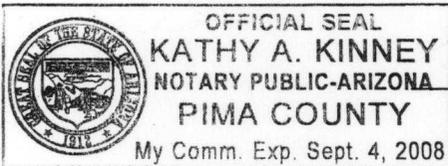

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

Ch
2.20.08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/14/08 TIME 9:51

INTOXILYZER SERIAL # 80-001199 LOCATION Team 2

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
02/14/2008 09:52:43

2 of 7 SC
Ch
2-20-08

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:53:30
Diagnostic Test	Pass	09:53:57
Air Blank	0.000	09:54:25
0.100 Cal Check	0.098	09:54:43
Air Blank	0.000	09:55:13
Subject Test	0.000	09:55:38
Air Blank	0.000	09:56:05
Five Minute Wait	ABT*	
Air Blank	0.000	09:56:46

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
02/14/2008 09:58:47

3 of 75c
M
2-22-08

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:59:36
Air Blank	IPS**	09:59:44

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 SC

OK
2-20-08

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
02/14/2008 10:01:51

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:02:41
Diagnostic Test	Pass	10:03:08
Air Blank	0.000	10:03:36
0.100 Cal Check	0.097	10:03:55
Air Blank	0.000	10:04:24
Subject Test	0.111*	10:07:32
Air Blank	0.000	10:08:04
Air Blank	0.000	10:08:31
Subject Test	XXX**	10:08:56
Air Blank	0.000	10:09:29

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence


INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
02/14/2008 10:09:53

5 of 75C
CU
2-20-08

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:10:20
Cal Check	0.096	10:10:39
Air Blank	0.000	10:11:08

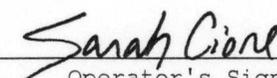

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

OK
2-20-08

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 2/14/08 TIME 9:51
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
02/14/2008 09:51:06

7 of 7 sc

ll
2-20-08

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:51:33
Cal Check	0.097	09:51:52
Air Blank	0.000	09:52:21

Sarah Cione

Operator's Signature