



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
2-19-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/13/08 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 2/14/08 and 3/13/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
 The instrument was returned to service on —.

Sarah Cione

Sarah L. Cione


2/19/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 19th day of February, 2009.

A.M. Hernandez

Notary Public

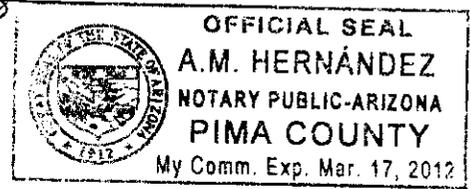

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7sc
TS4
3.17.08
TS4 2.14.09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/13/08 TIME 10:09

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

- (✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 sc
- (✓) 2. Error recognition logic system functioning. 100 sc AC.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system. 100 sc AC.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

On 2/19/09, I changed the alcohol-free subject test result from 0.100 to 0.000, as printed on page 2 of 7. This was done to correct a clerical error when I wrote the result of 0.100 on 3/13/08 sc

SIGNATURE

Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
03/13/2008

10:11:33

2 of 7sc
TS4

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:12:24
Diagnostic Test	Pass	10:12:52
Air Blank	0.000	10:13:19
0.100 Cal Check	0.098	10:13:38
Air Blank	0.000	10:14:07
Subject Test	0.000	10:14:24
Air Blank	0.000	10:14:52
Five Minute Wait	ABT*	
Air Blank	0.000	10:15:32

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
03/13/2008

10:17:39

3 of 7sc
TBY

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:18:26
Diagnostic Test	Pass	10:18:53
Air Blank	RFI*	10:19:02
Air Blank	IPS**	10:19:10

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
03/13/2008

10:19:27

4 of 7sc
134

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:20:12
Diagnostic Test	Pass	10:20:39
Air Blank	0.000	10:21:07
0.100 Cal Check	0.097	10:21:26
Air Blank	0.000	10:21:55
Subject Test	0.034*	10:25:03
Air Blank	0.000	10:25:36
Air Blank	0.000	10:26:03
Subject Test	XXX**	10:26:29
Air Blank	0.000	10:27:02

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
03/13/2008

10:27:49

5 of 75c
PSY

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:28:16
Cal Check	0.097	10:28:35
Air Blank	0.000	10:29:04

Sarah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 75c
TB4
3.17.2008

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/13/08 TIME 10:09
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
03/13/2008

10:09:42

Standard Lot#: 719801I
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:10:09
Cal Check	0.097	10:10:28
Air Blank	0.000	10:10:57

Sarah Cione

Operator's Signature

7 of 7sc

D4