



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

TJH  
 5-11-2009

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/7/09 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 4/9/09 and 5/7/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


5/7/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 7th, day of May, 2009.

A.M. Hernandez

Notary Public

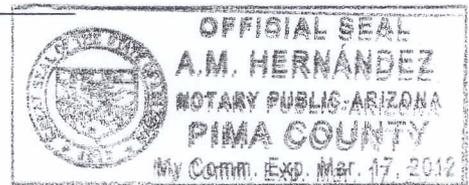

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
5.11.2009

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/7/09 TIME 10:26

INTOXILYZER SERIAL # 80-001163 LOCATION Vehicle # 868 - C15

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7 sc

154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
05/07/2009 10:28:00

-----  
Standard Lot#: 8043021  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

| Test             | g/210L | Time     |
|------------------|--------|----------|
| Air Blank        | 0.000  | 10:28:55 |
| Diagnostic Test  | Pass   | 10:29:24 |
| Air Blank        | 0.000  | 10:29:51 |
| 0.100 Cal Check  | 0.100  | 10:30:10 |
| Air Blank        | 0.000  | 10:30:39 |
| Subject Test     | 0.000  | 10:30:54 |
| Air Blank        | 0.000  | 10:31:21 |
| Five Minute Wait | ABT*   |          |
| Air Blank        | 0.000  | 10:31:56 |

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 7 sc

134

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
05/07/2009 10:32:02

-----  
Standard Lot#: 8043021  
Last Changed By: C. CONLEY #34517  
-----

OAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | RFI*   | 10:32:29 |
| Air Blank | IPS**  | 10:32:39 |

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4 of 7<sub>sc</sub>

154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
05/07/2009 10:32:45

-----  
Standard Lot#: 8043021  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 10:33:30 |
| Diagnostic Test | Pass   | 10:33:58 |
| Air Blank       | 0.000  | 10:34:26 |
| 0.100 Cal Check | 0.100  | 10:34:44 |
| Air Blank       | 0.000  | 10:35:14 |
| Subject Test    | 0.088* | 10:38:22 |
| Air Blank       | 0.000  | 10:38:54 |
| Air Blank       | 0.000  | 10:39:22 |
| Subject Test    | XXX**  | 10:39:42 |
| Air Blank       | 0.000  | 10:40:15 |

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

5 of 7 sc

134

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
05/07/2009 10:40:31

-----  
Standard Lot#: 8043021  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 10:40:59 |
| Cal Check | 0.099  | 10:41:17 |
| Air Blank | 0.000  | 10:41:47 |

*Sarah Cione*  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

BH  
5-11-2009

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 5/7/09 TIME 10:26  
Intoxilyzer Serial # 80-001163 Location Vehicle # 868 - C15

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 101 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

7 of 7sc

1/24

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
05/07/2009 10:26:37

-----  
Standard Lot#: 8043021  
Last Changed By: C. CONLEY #34517

-----  
QAS: S. CIONE #50616  
TUCSON PD

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 10:27:04 |
| Cal Check | 0.101  | 10:27:23 |
| Air Blank | 0.000  | 10:27:52 |

-----  
*Sarah Cione*  
Operator's Signature