



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TSY
6-8-2009

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/4/09 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 5/7/09 and 6/4/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


6/4/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 4th day of June, 2009.

A.M. Hernandez

Notary Public

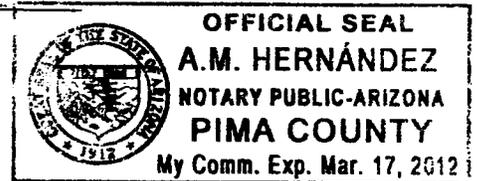

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TJH
6.8.2009

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 6/4/09 TIME 9:16
INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
06/04/2009 09:18:57

BM

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:19:42
Diagnostic Test	Pass	09:20:09
Air Blank	0.000	09:20:37
0.100 Cal Check	0.101	09:20:55
Air Blank	0.000	09:21:25
Subject Test	0.000	09:21:44
Air Blank	0.000	09:22:11
Five Minute Wait	ABT*	
Air Blank	0.000	09:22:45

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
06/04/2009

09:22:52

134

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:23:24
Air Blank	IPS**	09:23:32

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
06/04/2009

09:23:46

TS4

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:24:29
Diagnostic Test	Pass	09:24:57
Air Blank	0.000	09:25:24
0.100 Cal Check	0.100	09:25:43
Air Blank	0.000	09:26:12
Subject Test	0.154*	09:29:21
Air Blank	0.000	09:29:53
Air Blank	0.000	09:30:20
Subject Test	XXX**	09:30:38
Air Blank	0.000	09:31:11

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
06/04/2009 09:31:35

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

104

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:32:02
Cal Check	0.100	09:32:21
Air Blank	0.000	09:32:50

Sarah Cione

Operator's Signature

6 of 7sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSC
6-8-2009

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/4/09 TIME 9:16
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

7 of 7sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
06/04/2009 09:16:59

TS4

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:17:26
Cal Check	0.100	09:17:45
Air Blank	0.000	09:18:14

Sarah Cione
Operator's Signature